



+ Evidence in focus

PICO[◇] sNPWT clinical compendium

PICO[◇]
Single Use Negative Pressure
Wound Therapy System

September 2023

Smith+Nephew

Key for icons



Closed incisions



Chronic/open/
non-healing wounds



Health economics



Link to full
Evidence in focus
summary



Link to open
access reference



Link to published
abstract

Evidence pyramid



Level 1



Level 2



Level 3



Level 4



Level 5

Abbreviations

ABPI	Ankle-brachial pressure index	LS	Least squares	SSC	Surgical site complication
ASA	American Society of Anaesthesiologists	MTG	Medical technologies guidance	SSI	Surgical site infection
ASEPSIS	A quantitative scoring system used to identify and classify SSIs	NICE	National Institute for Health and Care Excellence	SWD	Surgical wound dehiscence
BIMA	Bilateral internal mammary artery	NNT	Number needed to treat	TAA	Total ankle arthroplasty
BHIS	Brompton and Harefield Infection Score	NPWT	Negative pressure wound therapy	TEWL	Transepidermal water loss
BMI	Body mass index	OR	Odds ratio	TKA	Total knee arthroplasty
CABG	Coronary artery bypass graft	POSAS	Patient Observer Scar Assessment Scale	tNPWT	Traditional negative pressure wound therapy
CDC	Centers for Disease Control and Prevention	PP	Per-protocol	VAS	Visual analogue scale
DFU	Diabetic foot ulcer	PU	Pressure ulcer	VLU	Venous leg ulcer
DSWI	Deep sternal wound infection	QALY	Quality adjusted life years	WUWHS	World Union of Wound Healing Societies
FEA	Finite element analysis	RCT	Randomised controlled trial		
ITT	Intention-to-treat	RRR	Relative risk reduction		
LoS	Length of stay	sNPWT	Single use negative pressure wound therapy		

Introduction

PICO[◇] Single Use Negative Pressure Wound Therapy (sNPWT) has a strong evidence base.

To date, **316**^{*} clinical publications (peer-reviewed manuscripts and conference abstracts) regarding PICO sNPWT have been identified (166 unique studies). This evidence compendium contains summaries of the most relevant publications; it does not include all the publications due to the volume of studies.

Levels of evidence*



60

RCTs, meta-analyses,
health economics
evaluations of RCTs

(+50 studies that note
PICO sNPWT studies)



19

Prospective
comparative
observational
studies



35

Retrospective
comparative
observational
studies



62

Case series
and case studies



90

Expert opinion,
case studies
or bench research

*Up to 18 May 2023. Evidence analysis report; EA/AWM/PICO/030/v3.


PICO^o sNPWT indication¹

PICO sNPWT is indicated for patients who would benefit from a suction device (NPWT) as it may promote wound healing via removal of low to moderate levels of exudate and infectious materials.

Appropriate wound types include:

- Closed surgical incisions
- Chronic
- Acute
- Traumatic
- Subacute and dehisced wounds
- Partial-thickness burns
- Ulcers (such as diabetic or pressure)
- Flaps and grafts


PICO sNPWT systems are suitable for use both in a hospital and homecare setting.

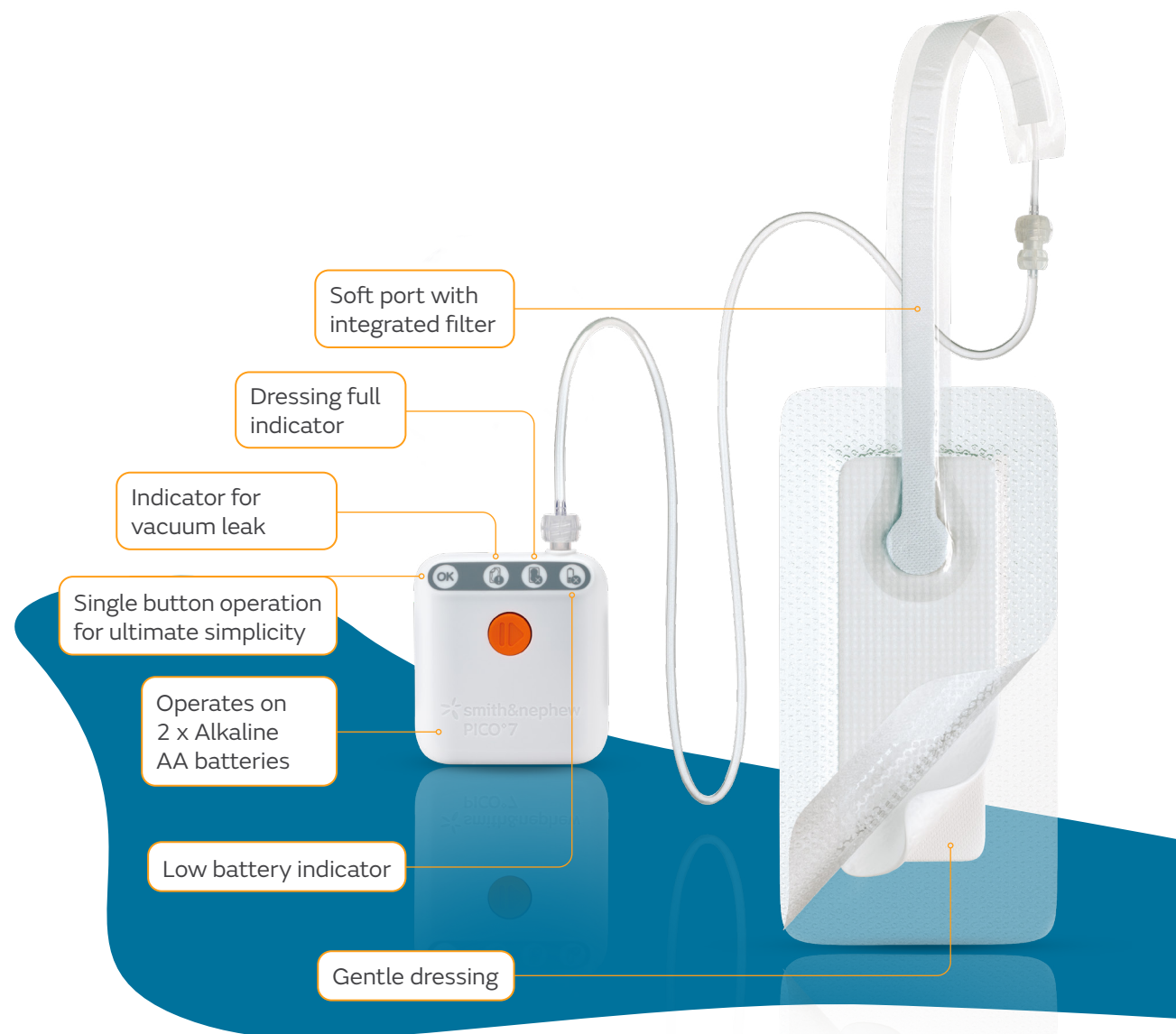

PICO[◇] sNPWT multilayer dressing with AIRLOCK[◇] Technology

Up to
80%
of the exudate
is lost
by evaporation⁹

Up to
20%
of the exudate
is absorbed
by dressing⁹

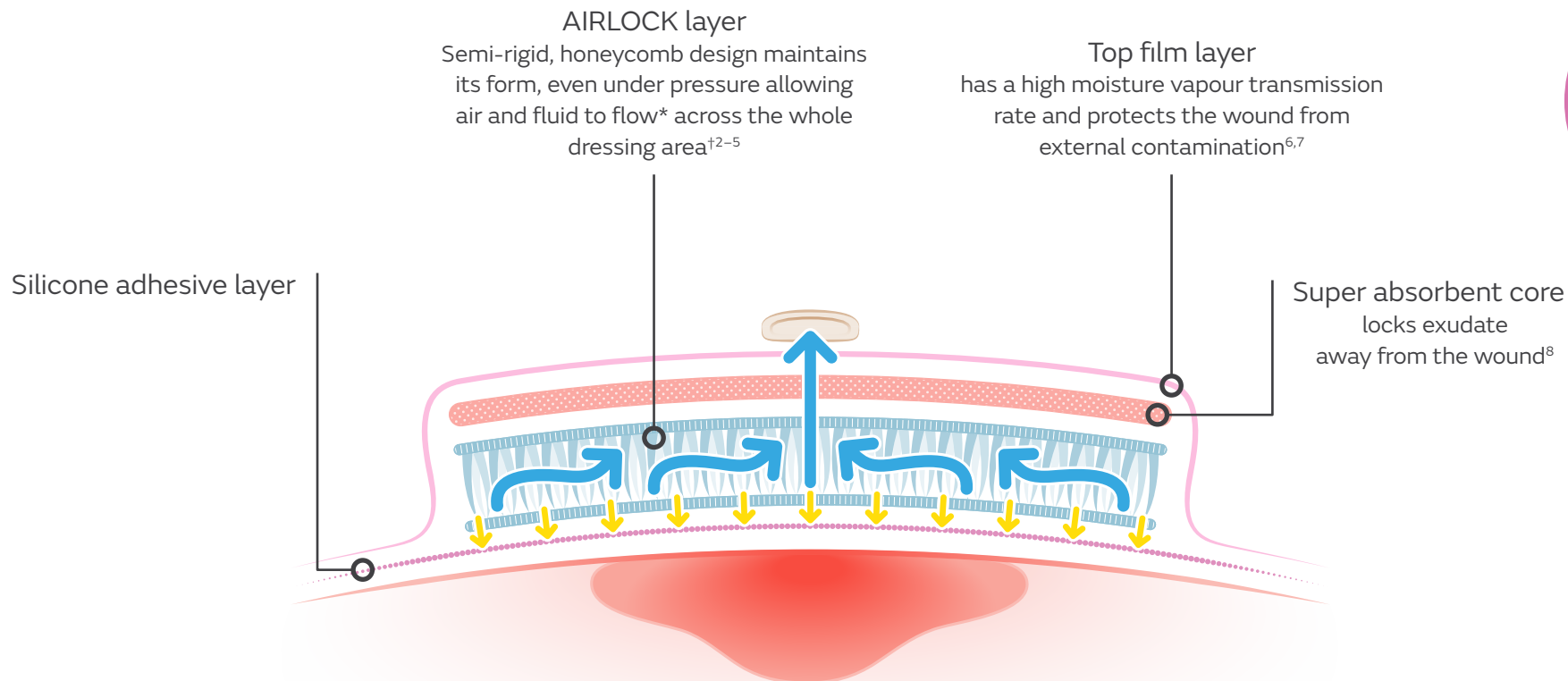

*As demonstrated in benchtop testing. †As demonstrated in ex vivo studies.

Consensus document/guidelines about prophylactic NPWT and PICO[◇] sNPWT for closed surgical incisions

World Health Organization

The World Health Organization recommends the use of prophylactic NPWT “in adult patients on primarily closed surgical incisions in high-risk wounds, for the purpose of the prevention of SSI, while taking resources into account.”¹⁰

World Union of Wound Healing Societies

WUWHS proposes NPWT is used in patients with closed surgical incisions who have intrinsic risk factors for SSCs or who have had a surgical procedure associated with higher incidence and/or higher consequence of SSCs.¹¹

The 2019 WUWHS Consensus Document on Wound Exudate: effective assessment and management, recognises the benefits of sNPWT in the management of closed surgical incisions:¹²

- Provides a barrier to external contamination^{11,12}
- Removes excess wound exudate¹²
- May aid healing by:^{11,12}
 - Reducing lateral tension across the closed incision
 - Improving lymphatic drainage
 - Reducing the risk of wound infection and separation (dehiscence)

National Institute for Health and Care Excellence (NICE)

NICE Medical technologies guidance: *PICO negative pressure wound dressings for closed surgical incisions* (MTG43).

NICE aims to improve health and social care in England through evidence-based guidance. NICE guidance helps people make efficient, cost-effective and consistent decisions about adopting new medical technologies. NICE guidance is internationally recognised.


NICE recommends that PICO sNPWT should be considered as an option for closed surgical incisions in patients who are at high risk of SSIs.¹³

In a review of data from 31 clinical studies (15 randomised controlled trials and 16 non-randomised comparative observational studies), NICE concluded that PICO sNPWT is associated with fewer SSIs and seromas compared with standard wound dressings. Cost modelling suggests that compared with standard wound dressings, PICO sNPWT provides extra clinical benefits at a similar overall cost with standard wound dressings.⁵

510k clearance

PICO sNPWT is the first system indicated in the US to aid in reducing the incidence of both deep and superficial incisional SSIs as well as post-operative seroma and dehiscence for high risk patients in Class I and Class II wounds.¹⁴

Introduction

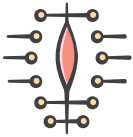

30% 
of wounds persist
for more than one year^{*35}

Unhealed wounds 
cost on average **135%**
more than wounds that have healed³⁶

PICO[®] sNPWT has been shown
to **significantly improve a variety
of healing-related outcomes in a range
of non-healing wounds indications**

NON-HEALING WOUNDS	
VLU, DFU	Kirsner R, et al. (2019)
	Kirsner RS, et al. (2020)
	Patel A, et al. (2019)
DFU	Sharpe A, et al. (2018)
VLU, PU	Hampton J. (2015)
Dehisced surgical wounds	Hughes J, et al. (2020)
Hard-to-heal wounds of various aetiologies	Hampton J, et al. (2022)
	Hurd T, et al. (2020)
	McCluskey P, et al. (2020)
	Dowsett C, et al. (2017)

 Click on the author
to navigate to study

^{*}Data representing the National Health Service in the UK.

1. Kirsner R, et al.



A prospective, randomised, controlled clinical trial on the efficacy of a single-use negative pressure wound therapy system, compared to traditional negative pressure wound therapy in the treatment of chronic ulcers of the lower extremities

Kirsner R, Dove C, Reyzelman A, Vayser D, Jaimes H. *Wound Repair Regen.* 2019;27(5):519–529.

Overview

- A randomised, controlled, multicentre study conducted at 16 centres in the USA and two centres in Canada to evaluate efficacy and safety of PICO[®] sNPWT or tNPWT to manage lower extremity ulcers (>4 weeks in duration)
- In total, 161 patients were included in the ITT population (101 VLUs; 60 DFUs) and were randomised to receive either PICO sNPWT (n=80) or tNPWT (n=81)
 - The PP population (non-inferiority analysis) included 115 patients (PICO sNPWT, n=64; tNPWT, n=51)

Results

- Reduction in wound area was significantly greater with PICO sNPWT than tNPWT in the PP population (88.7 vs 58.6% mean reduction; $p=0.003$) and the ITT population ($p<0.001$; Figure)
 - Significant LS mean reductions in wound area were also achieved with PICO sNPWT versus tNPWT in VLU (36.2%; $p=0.007$) and DFU (38.8%; $p=0.031$) subgroups
- Reductions in wound depth and volume in the PP and ITT populations (Figure) were also significantly greater with PICO sNPWT versus tNPWT ($p<0.02$, all comparisons)
- More patients had complete wound closure at 12 weeks with PICO sNPWT than with tNPWT (45 vs 22%; $p=0.002$; ITT population)
- Overall satisfaction with PICO sNPWT was significantly greater than with tNPWT

Conclusions

In patients with VLUs and DFUs, PICO sNPWT significantly reduced wound area, depth and volume compared with tNPWT; complete closure of lower extremity ulcers at 12 weeks was more frequent with PICO sNPWT than with tNPWT.

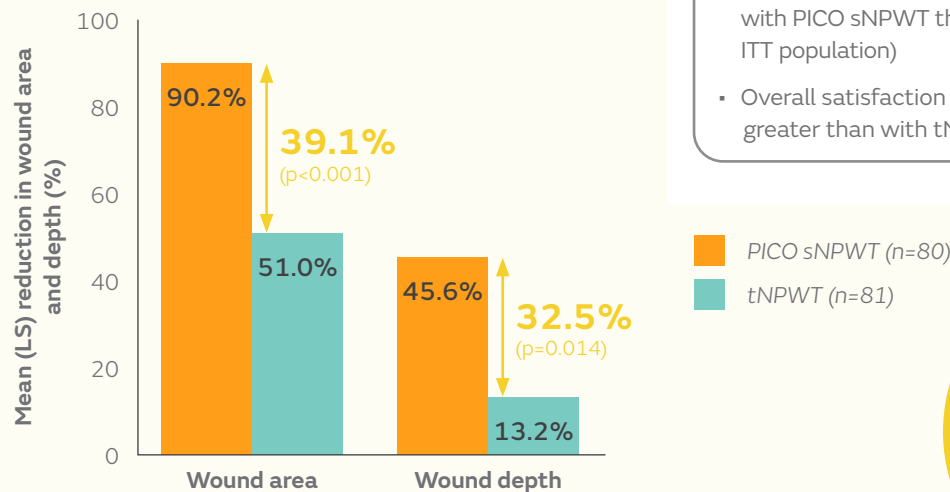


Figure. Percentage reductions from baseline in wound area and depth with PICO sNPWT and tNPWT at 12 weeks (ITT population; LS mean values)

51%
relative increase
in patients with complete
wound closure
at 12 weeks with PICO sNPWT
versus tNPWT ($p=0.002$)

2. Kirsner RS, et al.



A cost-effectiveness analysis comparing single-use and traditional negative pressure wound therapy to treat chronic venous and diabetic foot ulcers

Kirsner RS, Delhougne G, Searle RJ. *Wound Manag Prev.* 2020;66(3):30–38.

Overview

- A cost-effectiveness evaluation of PICO[®] sNPWT and tNPWT in treating lower extremity ulcers (US payer perspective)
 - Time horizons of 12 and 26 weeks were used to show the effect on wound closure
- Analysis of data from Kirsner, et al., 2019 and US National 2016 Medicare claims

Results

- For both ulcer types combined, switching from tNPWT to PICO sNPWT resulted in an estimated:
 - Expected cost saving per patient of \$7,756 at week 12 and \$15,749 at week 26
 - Decrease in total expected open ulcer weeks of 1.67 at week 12 and 5.31 at week 26
 - Increase in percentage of expected closed ulcers of 22.6% at week 12 and 31.0% at week 26
- Similar results were observed for VLU and DFUs when analysed separately

Conclusions

PICO sNPWT was estimated to be highly cost saving and reduced expected weeks to ulcer closure compared with tNPWT in patients with VLUs and DFUs, when analysed from a US payer perspective.

\$15,749 estimated

cost saving per patient with PICO sNPWT versus tNPWT at week 26



3. Patel A, et al.



Comparison of wound closure in chronic lower extremity ulcers between single use negative pressure wound therapy and traditional negative pressure wound therapy: a real-world analysis

Patel A, Delhougne G, Nherera L. Poster presented at: Wild on Wounds National Conference. September 11–14, 2019. Las Vegas, NV, USA.

Overview

- Retrospective cohort study to assess wound closure rates with PICO[®] sNPWT and tNPWT in a real-world setting in patients with DFUs and VLUs
 - PICO sNPWT: DFUs (n=84); VLUs (n=62)
 - tNPWT: DFUs (n=86); VLUs (n=60)

Results

- Compared with tNPWT, wound closure rates with PICO sNPWT were greater for all lower extremity ulcers (46.6 vs 34.9%; $p=0.043$)
 - Rates were also greater for DFUs and VLUs when analysed alone
- Compared with tNPWT, wounds treated with PICO sNPWT were 89% more likely to achieve closure ($p=0.042$)

Conclusions

Lower extremity ulcers (DFUs and VLUs) of patients treated with PICO sNPWT were more likely to achieve wound closure than those treated with tNPWT in this retrospective analysis of real-world outpatient wound clinic data.

4. Sharpe A, et al.



Using single use negative pressure wound therapy for patients with complicated diabetic foot ulcers: an economic perspective

Sharpe A, Myers D, Searle R. *Wounds UK*. 2018;14(4):89–93.

Overview

- UK case series of four patients using PICO 7 sNPWT to help manage complicated DFUs
- Patients and carers self-assessed the dressing status using the dressing-full indicator
 - PICO 7 sNPWT (n=4)

Results

- All four DFUs improved (mean ulcer area reduction, 49%), exudate levels were managed effectively and the frequency of dressing changes was reduced
- Total combined weekly clinician time saving using PICO 7 sNPWT was 279min (4hr 39min) for four patients
- Use of PICO sNPWT was estimated to release 13.5 clinician hours per patient on average over 12 weeks

Conclusions

Frequency of clinician visits and dressing changes were reduced by using PICO 7 sNPWT to help manage DFUs, improving service delivery with potential efficiency savings compared with prior practice.

5. Hampton J, et al.



Providing cost-effective treatment of hard-to-heal wounds in the community through use of NPWT

Hampton J. *Br J Community Nurs*. 2015;S14 (Suppl Community Wound Care): S16–S20.

Overview

- Cohort case study involving patients with hard-to-heal VLU and PUs treated in the community setting for >6 weeks
- Patients received PICO[®] sNPWT for 2 weeks followed by standard treatment appropriate for each wound
 - PICO sNPWT (n=9)

Results

- Average weekly reduction in wound size was 21%
- With PICO sNPWT target wound size was achieved on average 10 weeks earlier than predicted with standard treatment
- In wounds that responded, wound size reduction was 6 times faster than predicted with standard treatment
- Mean savings of DKK 6,670 (€895)* per patient using PICO sNPWT compared with prior standard treatment

Conclusions

Use of PICO sNPWT for 2 weeks helped to kick start the healing of chronic hard-to-heal wounds, which resulted in faster overall rates of healing and reduced costs compared with previous standard treatment.

*Exchange rate 1 EUR = 7.45550 DKK as of May 19, 2020.

6. Hughes J, et al.



The burden of dehisced wounds in the community: using early results from a multi-centre service evaluation to propose a standard of care to improve patient outcomes and safeguard woundcare budgets

Hughes J, Costello M, Belshaw M, Horton H, Styche T. *Br J Health Care Manag.* 2020;27:16–25.

Overview

- Analysis of a subset of 34 dehisced surgical wounds from a service evaluation of PICO[®] sNPWT
- The service evaluation adopted a previously published pathway for hard-to-heal wounds

Results

- After implementation of the PICO sNPWT pathway, 18 of 34 wounds (53%) healed within 12 weeks
 - Mean time to healing was 6.1 weeks
- Mean dressing change frequency reduced with use of PICO sNPWT from 4.7 to 3.2 times per week and remained at 3.3 times per week after returning to standard care
- Estimated savings of £16,577 for total wound care treatment over 12 weeks with PICO sNPWT versus standard care
 - Nursing time was reduced by 513 hours using PICO sNPWT compared with standard care

Conclusions

Use of PICO sNPWT as part of a pathway for hard-to-heal wounds helped to support healing of dehisced surgical wounds, as well as reduce estimated total costs and release nursing time compared with prior standard care in this service evaluation.



53%

of dehisced wounds healed within 12 weeks



513

nursing hours saved with PICO sNPWT versus standard care



21.6%

estimated cost reduction with PICO sNPWT use versus standard care (from £76,828 to 60,251)

7. Hampton J, et al.



Multi-centre, international practice-based evidence using PICO[®] single-use negative pressure wound therapy: challenging current behaviours in wound care practice

Hampton J, Meagher H, Sharpe A, et al. *Wounds International*. 2022;13(2):46–53.

Overview


- In service evaluation involving patients with non-healing wounds of various aetiologies:
 - Wounds were predominantly static or had minimal progression towards healing
 - Mean wound duration was 26.5 weeks
- 323 patients were treated with PICO sNPWT as part of a predetermined clinical pathway
- An economic model compared outcomes with the likely outcomes had PICO sNPWT not been employed

Results

- Within 12 weeks of initiation of PICO sNPWT:
 - 52% of the wounds healed
 - Dressing change frequency reduced by a third (vs frequency before PICO sNPWT; 3.0 vs 4.7 times per week)
 - Costs were estimated to have reduced by 30% to £651 per patient

Conclusions

With use of PICO sNPWT in a pathway for non-healing wounds, more than half of wounds had healed in 12 weeks and dressing change frequency was reduced compared with standard care. Overall wound care costs were estimated to reduce by around 30%.

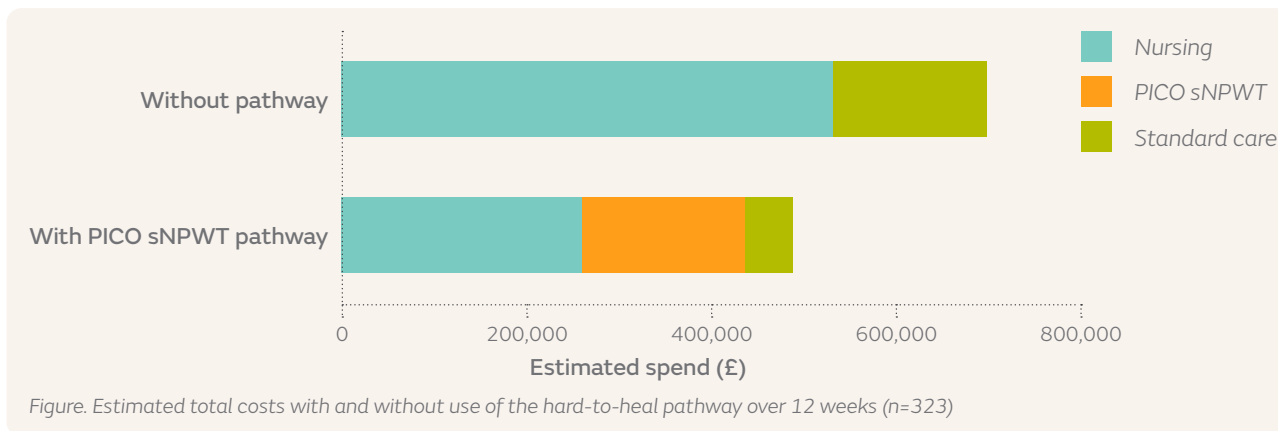

8. Hurd T, et al.



Single use negative pressure wound therapy (sNPWT) in the community management of chronic open wounds deeper than 2cm

Hurd T, Gilchrist B. Poster presented at: Symposium on Advanced Wound Care/WHS Annual Meeting. July 24–26, 2020; virtual conference.

Overview

- Retrospective two-year analysis of the healing of chronic open wounds >2cm deep (DFUs, VLU, PUs and dehiscent surgical wounds) in the home or community care setting following introduction of an integrated care bundle including PICO[®] sNPWT compared with standard care
 - PICO sNPWT (409 wounds; patients were significantly older with higher comorbidity score, both $p < 0.001$)
 - Standard care (2,242 wounds)

Results

- Use of PICO sNPWT to manage chronic open wounds >2cm compared with standard care resulted in:
 - Shorter mean healing times (46% relative reduction; 11.5 days)
 - Longer mean time between dressing changes (3.23 days)

Conclusions

PICO sNPWT may help reduce healing times and frequency of dressing changes in chronic open wounds >2cm deep compared with standard care.

9. McCluskey P, et al.



McCluskey P, Brennan K, Mullan J, et al. JCN. 2020;34:36–43.

Overview

- A service evaluation at seven centres in Northern Ireland and the Republic of Ireland
- Wound healing and health economic impact (in UK sterling and Euros) of using PICO[®] sNPWT versus standard care on hard-to-heal wounds over 12 weeks (or until healing) were assessed
- Median wound duration was 3–6 months; 36 wounds were included
- Eligible patients had:
 - Wounds >6 weeks in duration with no signs of clinical infection
 - <10% per week wound area reduction over 4 weeks
 - No NPWT in the last 6 weeks or contraindications for NPWT
 - ABPI >0.8 and <1.3 for VLUs


Results

- Using PICO sNPWT, 20 of 36 wounds healed within 12 weeks (55.6%)
 - Mean healing time was 6.95 weeks
- Wound healing rate was greater for wounds with <3 months duration than those with ≥3 months duration (84.6 vs 71.4%; $p=0.0125$; Figure)
- Improvements in mean wound area per week with PICO sNPWT (-16.8%) continued after use (-18.9%)
- Dressing changes per week were less frequent with PICO sNPWT versus standard care (1.75 vs 3.56 changes; $p<0.001$)
 - They were also less frequent in the post PICO sNPWT phase (1.95 vs 3.56 changes per week; $p<0.001$)
- Use of PICO sNPWT was predicted to reduce costs versus standard care (Figure):
 - Total costs by 25% (£15,467) and 21% (€12,001)
 - Nursing resource costs by 59% (£31,494 and €27,517)

Impact of a single-use negative pressure wound therapy system on healing

Conclusions

In patients with hard-to-heal wounds, PICO sNPWT was most effective for wounds of <3 months in duration. It helped to reduce dressing change frequency and was predicted to reduce nursing resource costs compared with standard care.

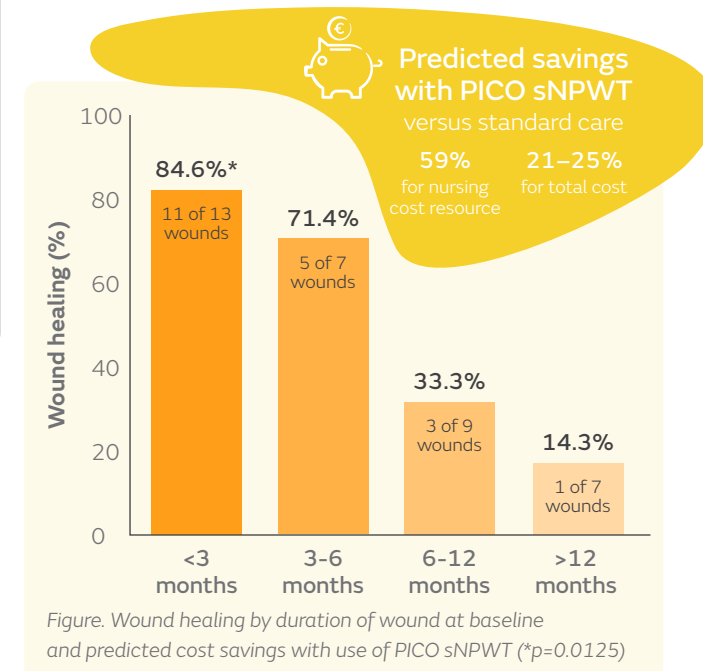

10. Dowsett C, et al.



Use of PICO[®] to improve clinical and economic outcomes in hard-to-heal wounds

Dowsett C, Hampton J, Myers D, Styche T. *Wounds International*. 2017;8(2):52–58.

Overview

- A prospective cohort study of 52 hard-to-heal wounds of varied aetiology and duration treated according to the PICO sNPWT pathway
 - Patients were switched from standard care to treatment with PICO sNPWT at week 0 for at least two weeks

Results

- During PICO sNPWT treatment, wound area reduced by 13.4% more per week than pre-PICO sNPWT ($p=0.006$)
- After the PICO sNPWT phase, wound area reduced by 9.6% more per week than pre-PICO sNPWT ($p=0.001$)
- PICO improved the trajectory of wounds of over 1 year, and healing rates were almost three times greater in wounds of <3 months duration (94.1 vs 33.3%)
- Implementing the PICO sNPWT pathway was estimated to reduce total costs by 33.1% (£50,000) and release 119 nursing days over 26 weeks compared with predictions for standard care

Conclusions

PICO sNPWT helped to significantly improve the healing trajectory of hard-to-heal wounds, resulting in reduced estimated costs and nursing time compared with previous standard care.




Estimated released
nursing days with PICO sNPWT
compared with predictions for standard care



33.1% estimated
cost reduction with PICO sNPWT
compared with predictions for standard care

Hard-to-heal pathway: when to use PICO[◇] sNPWT

Figure adapted from Dowsett C, et al. (2017)

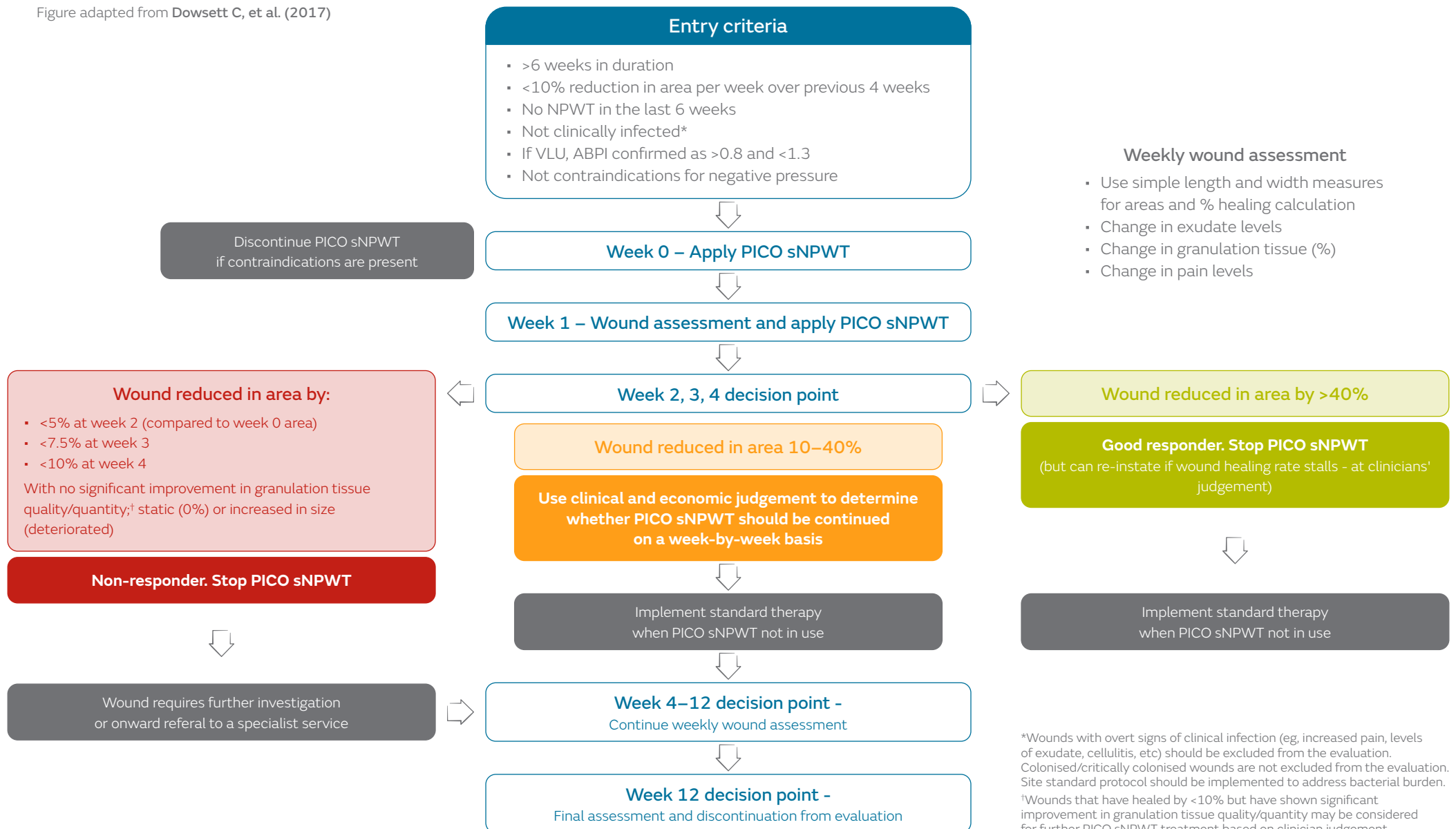

*Wounds with overt signs of clinical infection (eg, increased pain, levels of exudate, cellulitis, etc) should be excluded from the evaluation. Colonised/critically colonised wounds are not excluded from the evaluation. Site standard protocol should be implemented to address bacterial burden.

†Wounds that have healed by <10% but have shown significant improvement in granulation tissue quality/quantity may be considered for further PICO sNPWT treatment based on clinician judgement.

PICO[◇] sNPWT provides NPWT at -80mmHg

Clinical guidelines and consensus groups recommend therapeutic negative pressure levels of **-50 to -150mmHg** for wound care^{37,39}

High or low negative pressure levels are advised depending on exudate levels, wound type and pain experienced by patients^{37,39,40}

PICO sNPWT consistently delivers negative pressure at **-80mmHg**, a level sufficient to manage most wounds with low to moderate exudate⁴⁰

‘There is seldom any reason to use a negative pressure greater than -80mmHg, but as the drainage of exudate may be improved at -125mmHg, this pressure level could be used during the initial treatment of highly exuding wounds.’⁴⁰


Click here
for more details

For detailed product information, including indications for use, contraindications, precautions and warnings, please consult the product's applicable Instructions for Use (IFU) prior to use.

Advanced Wound Management
Smith & Nephew Medical Ltd
101 Hessle Road
Hull HU3 2BN, UK

www.smith-nephew.com
T +44 (0) 1482 225181
F +44 (0) 1482 328326

◇Trademark of Smith+Nephew
All Trademarks acknowledged
©September 2023 Smith+Nephew
AWM-AWD-26268 | GMC1149b

Developed by Evidence Communications,
Global Clinical & Medical Affairs

www.smith-nephew.com

1. PICO® Single Use Negative Wound Pressure. Available at: <https://www.smith-nephew.com/new-zealand/healthcare/products/product-types/negative-pressure-wound-therapy/pico/>. Accessed June 11, 2020. 2. Casey C. Consistent delivery of therapeutic negative pressure levels by a single use negative pressure wound therapy system (SNPWT)* in a wound model. Paper presented at: EWMA; 2019; Gothenburg, Sweden. 3. Smith+Nephew 2019. Internal Report. RD/19/006. 4. Smith+Nephew 2019. Report. DS/19/211/R. 5. Smith+Nephew 2021. Internal Report. DS/19/211/R - Part B. 6. Smith+Nephew 2018. Internal Report. DS/18/153/R. 7. Smith+Nephew 2020. Internal Report. 2001002. 8. Smith+Nephew 2018. Internal Report. DS.18.260.R. 9. Malmjö M, Huddleston E, Martin R. Biological effects of a disposable, canisterless negative pressure wound therapy system. *ePlasty*. 2014;14:e15. 10. WHO Global Guidelines for the Prevention of Surgical Site Infection. 2016. Available at: <https://apps.who.int/iris/handle/10665/277399>. Accessed August 30, 2023. 11. World Union of Wound Healing Societies (WUWHS) Consensus Document. Closed surgical incision management: understanding the role of NPWT. *Wounds International*, 2016. Available at: <https://woundsinternational.com/world-union-resources/closed-surgical-incision-management-understanding-the-role-of-npwt/>. Accessed August 30, 2023. 12. World Union of Wound Healing Societies (WUWHS) Consensus Document. Wound exudate: effective assessment and management. *Wounds International*, 2019. Available at: <https://woundsinternational.com/world-union-resources/wuwhs-consensus-document-wound-exudate-effective-assessment-and-management/>. Accessed August 30, 2023. 13. NICE Medical technologies guidance: PICO sNPWT negative pressure wound therapy for closed surgical incision wounds (MGT43). Available at: <https://www.nice.org.uk/guidance/mtg43>. Accessed June 11, 2020. 14. Smith+Nephew's PICO 7 and PICO 14 Negative Pressure Wound Therapy Systems are the first systems indicated to aid in reducing the incidence of both deep and superficial incisional surgical site infections and dehiscence. Available at: <http://www.smith-nephew.com/en-gb/news/2022/01/24/20220124-pico-7-and-pico-14-negative-pressure-wound-therapy-systems>. Accessed September 8, 2023. 15. Tanner J, Khan D, Aplin C, Ball J, Thomas M, Bankart J. Post-discharge surveillance to identify colorectal surgical site infection rates and related costs. *J. Hosp. Infect.* 2009;72(3):243–50. 16. Petrosillo N, Drapeau CM, Nicastrì E, Martini L, Ippolito G, Moro ML. Surgical site infections in Italian hospitals: a prospective multicenter study. *BMC Infectious Diseases*. 2008;8(1):34. 17. Wick EC, Gibbs L, Indorf LA, Varma MG, García-Aguilar J. Implementation of quality measures to reduce surgical site infection in colorectal patients. *Dis. Colon Rectum*. 2008;51(7):1004–1009. 18. Smith RL, Bohl JK, McElearney ST, et al. Wound infection after elective colorectal resection. *Ann. Surg.* 2004;239(5):599. 19. Yoshida J, Shinohara M, Ishikawa M, Matsuo K. Surgical site infection in general and thoracic surgery: surveillance of 2,663 cases in a Japanese teaching hospital. *Surg. Today*. 2006;36(2):114–8. 20. Serra-Aracil X, Espin-Basany E, Biondo S, Guirao X, Orrego C, Sitges-Serra A. Surgical site infection in elective operations for colorectal cancer after the application of preventive measures. *Arch. Surg.* 2011;146(5):606–12. 21. Olsen MA, Chu-Ongsakul S, Brandt KE, Dietz JR, Mayfield J, Fraser VJ. Hospital-associated costs due to surgical site infection after breast surgery. *Arch. Surg.* 2008;143(1):53–60. 22. Irwin GW, Boundouki G, Fakim B, et al. Negative Pressure Wound Therapy reduces wound breakdown and implant loss in prepectoral breast reconstruction. *Plast Reconstr Surg Glob Open*. 2020;8:e2667. 23. McIntosh J, O'Donoghue JM. Therapeutic mammoplasty – a systematic review of the evidence. *Eur. J. Surg. Oncol.* 2012;38(3):196–202. 24. Vincent SM, Gallagher M, Johnston A, Djohan R, Varzalis M, Sugrue M. The keys to optimising breast wounds: a meta-analysis. *Advances in Breast Cancer Research*. 2019;8(3):87–111. 25. Nelson JA, Fischer JP, Chung CU, et al. Obesity and early complications following reduction mammoplasty: an analysis of 4,545 patients from the 2005–2011 NSQIP datasets. *J. Plast. Surg. Hand Surg.* 2014;48(5):334–339. 26. Jenks PJ, Laurent M, McQuarry S, Watkins R. Clinical and economic burden of surgical site infection (SSI) and predicted financial consequences of elimination of SSI from an English hospital. *J. Hosp. Infect.* 2014;86(1):24–33. 27. Turtiainen J, Hakala T, Hakkarainen T, Karhukorpi J. The impact of surgical wound bacterial colonization on the incidence of surgical site infection after lower limb vascular surgery: a prospective observational study. *Eur. J. Vasc. Endovasc. Surg.* 2014;47(4):411–7. 28. Totty JP, Moss JW, Barker E, et al. The impact of surgical site infection on hospitalisation, treatment costs, and health-related quality of life after vascular surgery. *Int. Wound J.* 2021;18(3):261–8. 29. Audu CO, Columbo JA, Sun SJ, et al. Variation in timing and type of groin wound complications highlights the need for uniform reporting standards. *J. Vasc. Surg.* 2019;69(2):532–43. 30. NHS Digital. Maternity Services Monthly Statistics May 2020, experimental statistics. Available at: <https://digital.nhs.uk/data-and-information/publications/statistical/maternity-services-monthly-statistics/may-2020>. Accessed July 30, 2023. 31. Temming LA, Raghuraman N, Carter EB, et al. Impact of evidence-based interventions on wound complications after cesarean delivery. *Am J Obstet Gynecol.* 2017;217(4):449.e1–449.e9. 32. Patel VP, Walsh M, Sehgal B, Preston C, DeWal H, Di Cesare PE. Factors associated with prolonged wound drainage after primary total hip and knee arthroplasty. *J. Bone Joint Surg. Am.* 2007;89(1):33–38. 33. Poultsides LA, Ma Y, Della Valle AG, Chiu YL, Sulco TP, Memtsoudis SG. In-hospital surgical site infections after primary hip and knee arthroplasty – incidence and risk factors. *J. Arthroplasty*. 2013;28(3):385–389. 34. Tucker A, Walls A, Leckey B, et al. Postdischarge unscheduled care burden after lower limb arthroplasty. *J. Arthroplasty*. 2018;33(9):2745–2751. 35. Guest JF, Fuller GW, Vowden P. Cohort study evaluating the burden of wounds to the UK's National Health Service in 2017/2018: update from 2012/2013. *BMJ Open*. 2020 Dec 22;10(12):e045253. 36. Guest JF, Ayoub N, McIlwraith T, Uchegbu I, Gerrish A, Weidlich D, Vowden K, Vowden P. Health economic burden that different wound types impose on the UK's National Health Service. *Int Wound J.* 2017 Apr;14(2):322–330. 37. Apelqvist J, Willy C, Fagerdahl AM, et al. Negative Pressure Wound Therapy – overview, challenges and perspectives. *J Wound Care* 2017; 26: 3, Suppl 3, S1–S113. 38. Sugrue M, Ciprandi G, Djohan R, et al. World Union of Wound Healing Societies (WUWHS) Consensus Document. Closed surgical incision management: Understanding the role of NPWT. Available at: [www.woundsinternational.com/wuwhs/view/consensus-document-closed-surgical-incision-management-understanding-the-role-of-npwt](https://woundsinternational.com/wuwhs/view/consensus-document-closed-surgical-incision-management-understanding-the-role-of-npwt). Accessed July 30, 2023. 39. Birke-Sorensen H, Malmjö M, Rome P, et al. Evidence-based recommendations for negative pressure wound therapy: treatment variables (pressure levels, wound filler and contact layer) – steps towards an international consensus. *J Plast Reconstr Aesthet Surg.* 2011;64 Suppl:S1–16. 40. Malmjö M, Borgquist O. NPWT settings and dressing choices made easy. *Wounds International*. 2010;1(3):1–6.