


+ Rediscover normal

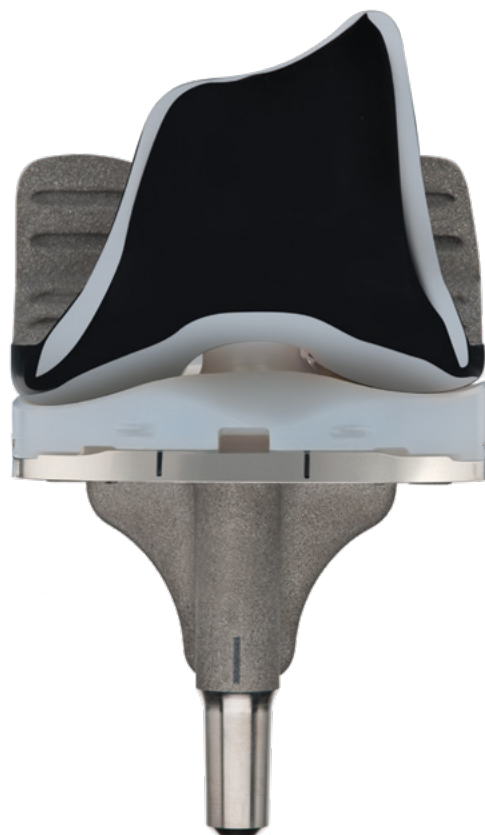

**Smith+Nephew**

JOURNEY<sup>◇</sup> II  
Total Knee Arthroplasty

**Product overview**

JOURNEY® II  
Total Knee Arthroplasty



# Rediscover normal: JOURNEY<sup>◇</sup> II TKA


JOURNEY II TKA is designed to restore normal shapes, position and motion<sup>1-5</sup> to help patients rediscover their normal through a smoother recovery<sup>\*\*6,7</sup>, improved function<sup>\*7-11</sup> and higher patient satisfaction.<sup>\*7-9,12</sup>

- + Normal shapes
- + Normal position
- + Normal motion




# Patient dissatisfaction after total knee replacement

Although total knee replacement has proved to be a successful treatment for improving pain and function, patients still report unmet levels of satisfaction, experience functional limitations and demonstrate increased difficulties doing daily activities compared to a healthy cohort. These deficits in patient satisfaction are even higher with more athletic and demanding activities.<sup>13,14</sup>

**>50%**  report **some degree of limitation** to their functional activities.<sup>13</sup>

## Patient difficulties in daily activities post-TKR<sup>13</sup>

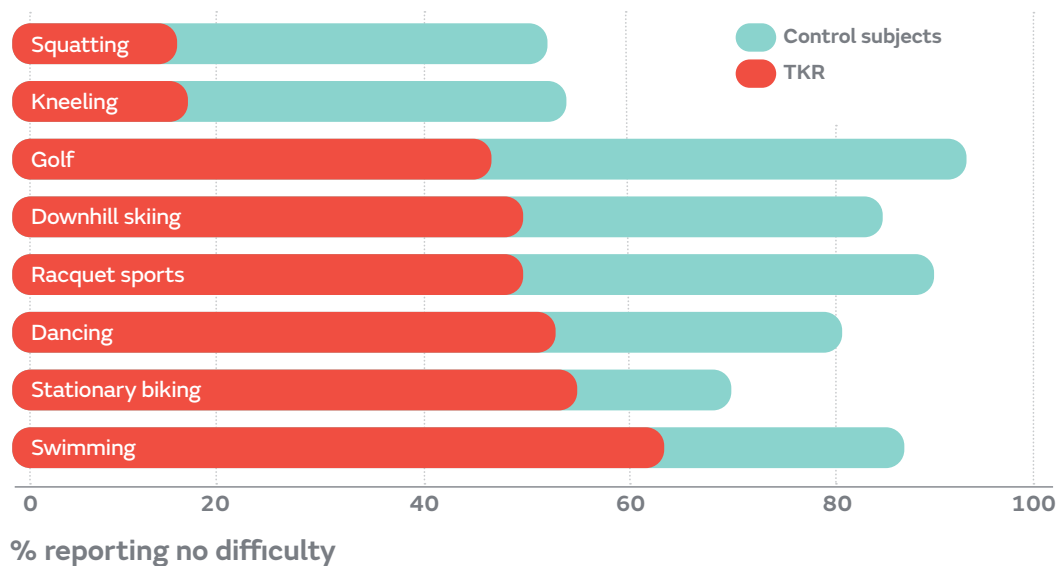

**20%**  of total knee replacement patients report **unmet levels of satisfaction**.<sup>14</sup>

# Implant design impacts satisfaction and performance

Through studying patient outcomes and performance, it has been observed that the motion, or kinematics, of the knee changes after total knee replacement. The differences found in the kinematics between normal knees and today's TKA designs have a direct impact on the feel, function, and performance of a patient's knee replacement.<sup>15-18</sup>

Today's TKA designs simply do not recreate the kinematics found in the normal knee, leading patients to experience functional limitations and feel their "artificial" knee.<sup>18</sup>



*The fact that many patients are unable to return to prior levels of function is partly attributable to differences in kinematics between the normal knee and current TKA designs<sup>17</sup>*


# 80%

of TKA patients report their knee feels **"Artificial"**<sup>18</sup>

# Normal knee kinematics

In order to understand patient dissatisfaction in total knees, it is important to understand what the normal knee looks like and how it works. The normal, healthy knee has unique shapes and positions that allow it to move and function the way it does.<sup>4,19</sup>

Studies on knee kinematics show that as the knee flexes, the femur externally rotates on the tibia. This external rotation is a result of the geometries of the medial and lateral sides of the tibia. The medial surface of the tibia is concave and restricts the amount of translation, causing the **medial femoral condyle to demonstrate a pivot-like motion**. However, the lateral surface of the tibia is convex and allows **the lateral femoral condyle to rollback posteriorly**.<sup>4,19</sup>

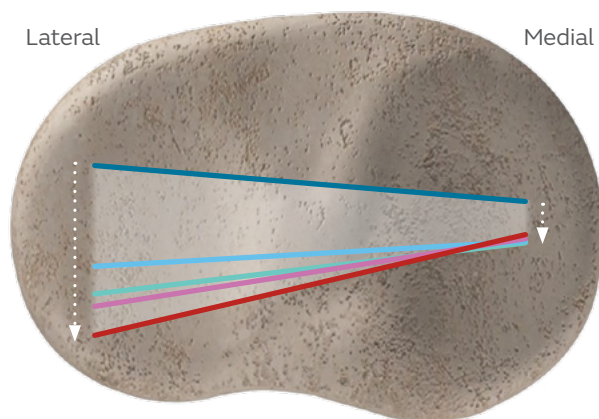


# Why do the kinematics change after total knee replacement?

## Implant design.

All TKA designs are essentially the same: symmetric and non-anatomic. These designs do not replicate the kinematics found in the normal knee because they do not replicate the shapes or positions of the normal knee.<sup>20-22</sup>

- Symmetric shapes and thicknesses (distal and posterior) of the femoral component lead to a non-anatomic, 0° joint line.<sup>20-22</sup>
- Femoral component sits in a non-anatomic posterior position on the tibia resulting in paradoxical motion and decreased muscular efficiency.<sup>20-22</sup>
- Lateral surface of the tibia insert has non-anatomical concave shape.<sup>20-22</sup>

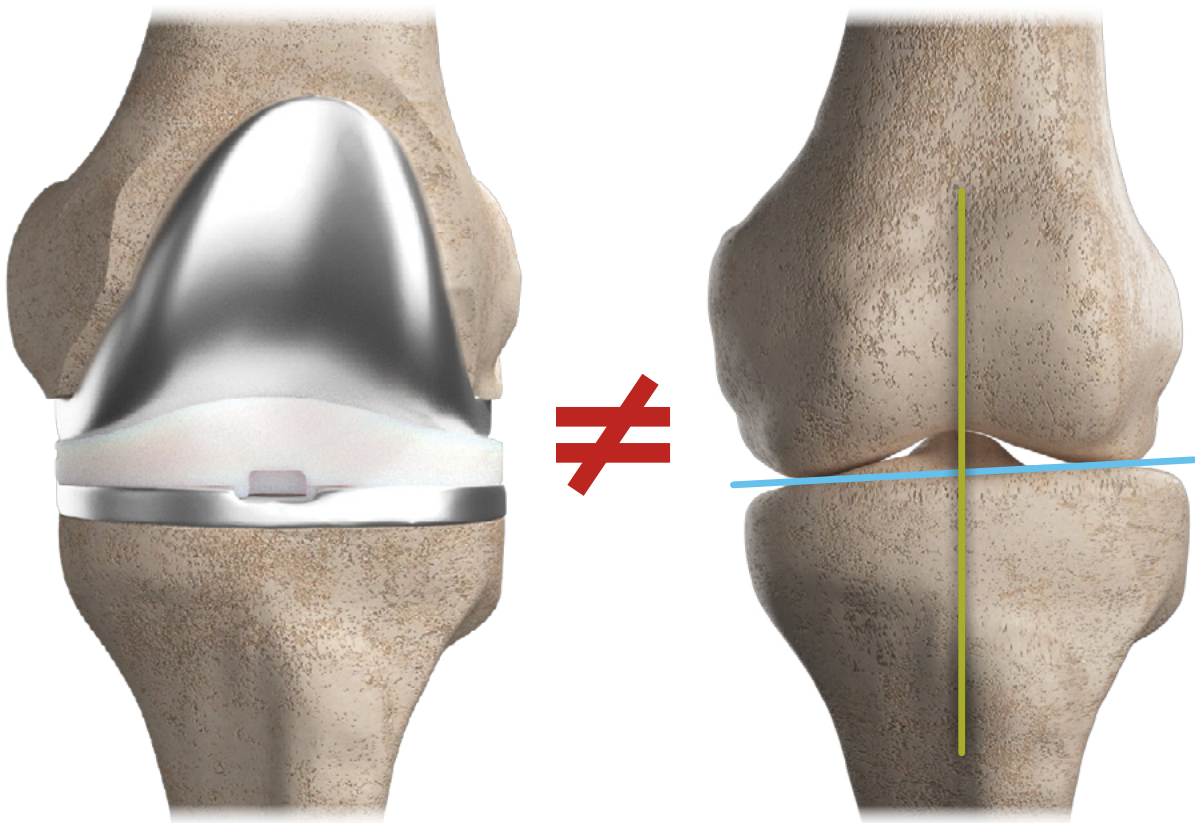


# Normal Knee vs Current TKA Shape

## Normal knee<sup>4,19</sup>

### Femur


The **medial condyle is thicker and rounder** than the lateral condyle

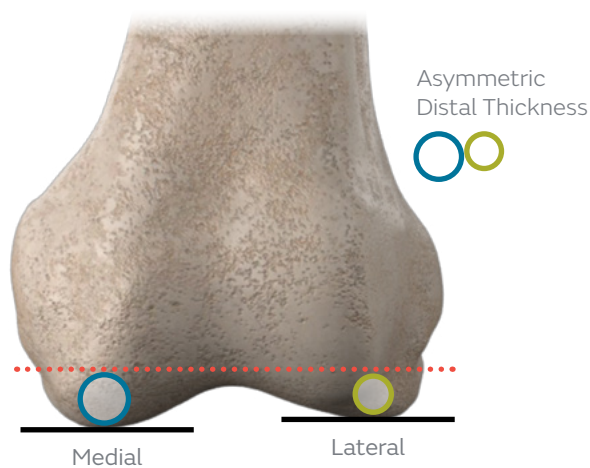

## Current TKA<sup>7,20-23</sup>


### Femur

The medial and lateral sides of the femoral component have **symmetric** distal and posterior condyles that are identical in thickness.

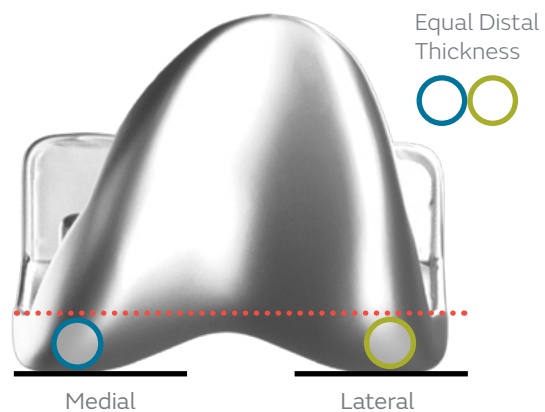


### Tibia

- Opposite of the femur, the lateral side of the tibia is thicker than the medial side
- Medial surface of the tibia is **concave** while the lateral surface is **convex**

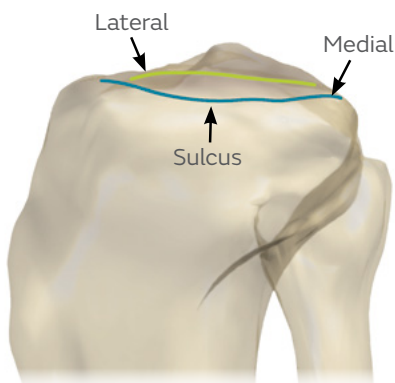

### Tibia

- Medial surface the **same thickness (symmetric)** as lateral surface
- Lateral surface of the tibia has **non-anatomical concave shape**

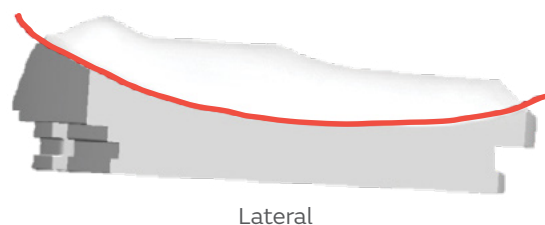




# Normal Knee vs Current TKA Position

## Normal knee<sup>4,19</sup>

### Joint line

When the asymmetric shapes of the femur and tibia are combined, they create a natural 3° varus joint line in the average knee.

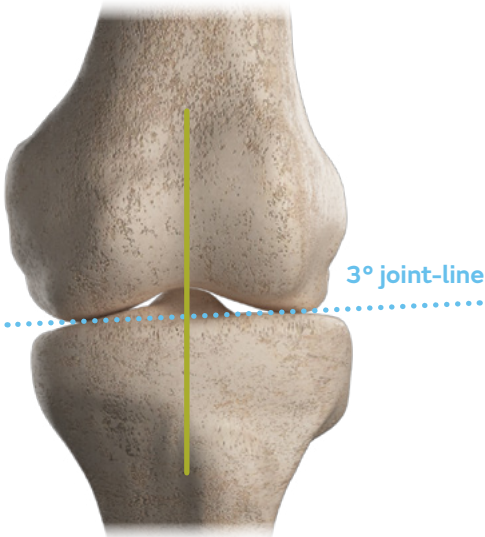

### A/P Position


The tibia has a **midline sulcus position** that causes the posterior femur to sit in-line with the posterior tibia with little to no overhang.



## Current TKA<sup>7,20-23</sup>


### Joint line

Non-anatomic, 0° joint line created by the symmetric shapes and thicknesses of the femur and tibia.

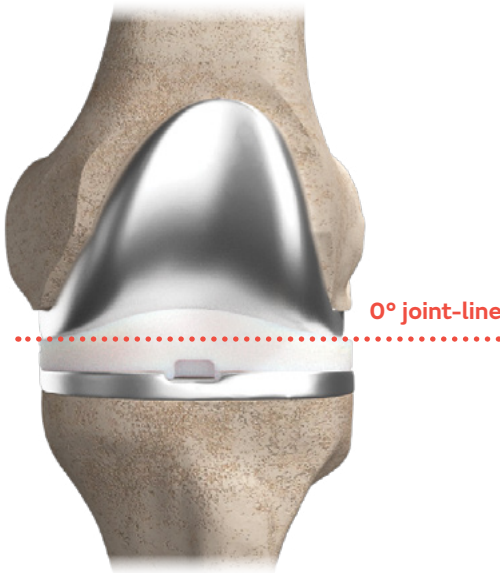


### A/P Position

Tibial insert has been designed with a more **posterior sulcus position** in order to get better flexion. This causes the femur to overhang the tibia.

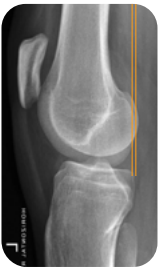

Normal

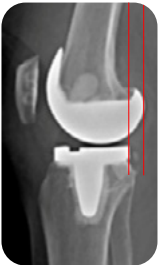

DePuy Attune™

# Normal Knee vs Current TKA Motion

Because TKA designs do not replicate the shapes and position of the normal knee, it is not possible to replicate the motion of the normal knee. These TKA designs experience:

- Paradoxical motion<sup>15,24</sup>
- Little to no external rotation throughout flexion<sup>15,24</sup>

## Normal knee<sup>4,19</sup>

### 0° (Full extension)

- Screw-home (5° femoral internal axial rotation)
- No posterior femoral overhang

### 1-90° (Mid flexion)

- Minimal medial translation (medial pivot)
- Lateral posterior rollback and external rotation

### 90-155° (Full flexion)

- Posterior femoral translation
- External axial rotation retained

## Current TKA<sup>20-23</sup>

### 0° (Full extension)



- No Screw-home to provide stability in extension
- Posterior sulcus position causes the femur to overhang the tibia posteriorly

### 1-90° (Mid flexion)

- **Paradoxical motion** is caused by the posterior position of the femoral component
- Femoral external axial rotation resisted by the concave surfaces of both the medial and lateral tibia (limited lateral rollback)


### 90-155° (Full flexion)

- Femoral component abnormally rotates internally and aligns with symmetric insert

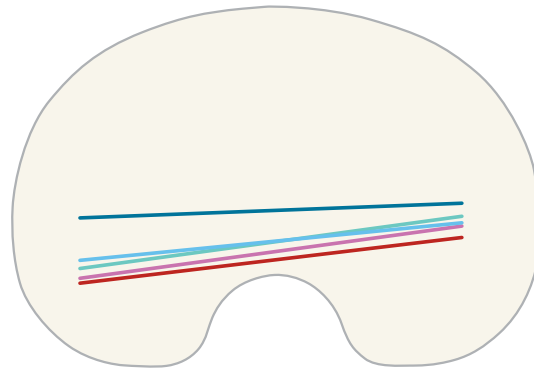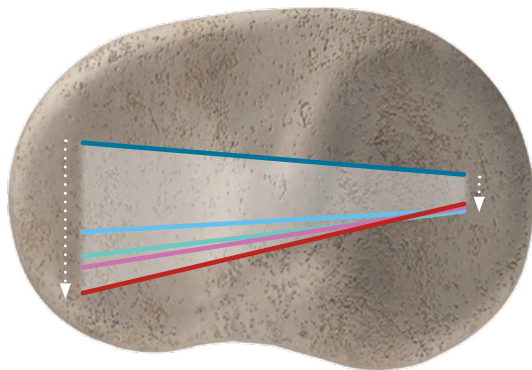


# Symmetric, non-anatomic shapes and positions lead to changes in:


## Gait cycle<sup>9,25-29</sup>

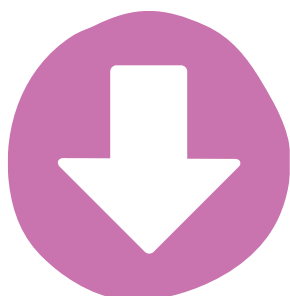

Velocity

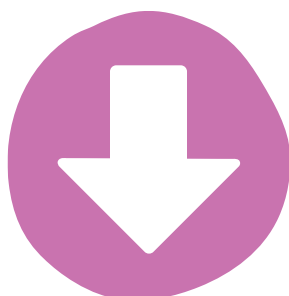

Stride length


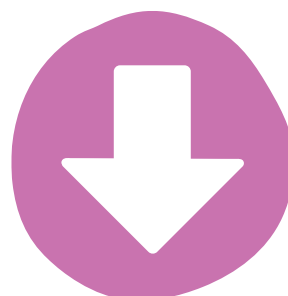


Max knee flexion  
during stance and  
swing phases




Mid-stance  
knee flexion  
(Quad avoidance gait)

## Muscle efficiency<sup>30</sup>

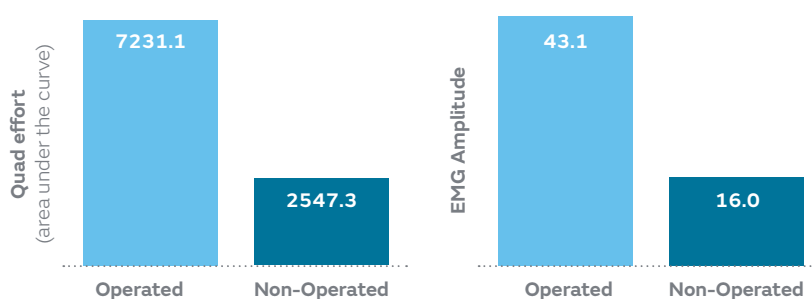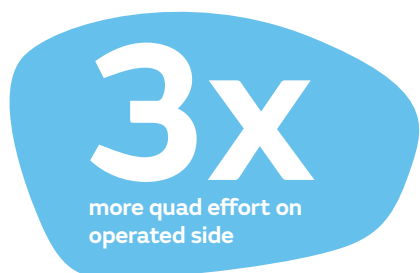

## Knee kinematics<sup>15,24</sup>

- More Posterior starting position than normal knee
- Paradoxical motion -> mid-flexion instability
- Little to no external rotation during flexion

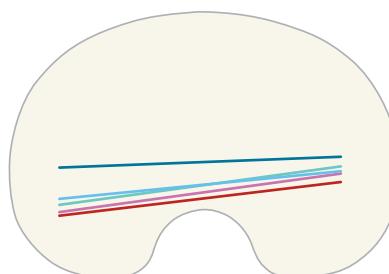

# Rediscover normal: JOURNEY<sup>◊</sup> II TKA

**Van Onsem et al.** stated, “Reproduction of optimal kinematic patterns during TKA could be instrumental in improving patient satisfaction after total knee replacement.”<sup>16</sup> The solution to providing patients with better overall satisfaction and functionality is to design an implant as close to the normal knee as possible. JOURNEY II was designed to do just that.

**Normal shapes + Normal position = Normal motion**


JOURNEY II TKA has been demonstrated to restore anatomical shape, position and motion.<sup>1-5</sup> This anatomical restoration can provide superior clinical outcomes and higher patient satisfaction.<sup>\*6-9,11,13</sup>



# Normal shapes<sup>4,5,8,31,32</sup>

**JOURNEY II TKA is designed to replicate the anatomic shapes found in the normal, healthy knee. These unique shapes include:**

- + Anatomic, asymmetric Femur/Tibia<sup>4,5,8,31,32</sup>
- + Concave **medial** tibial surface<sup>4,5,31</sup>
- + Convex **lateral** tibial surface<sup>4,5,31</sup>

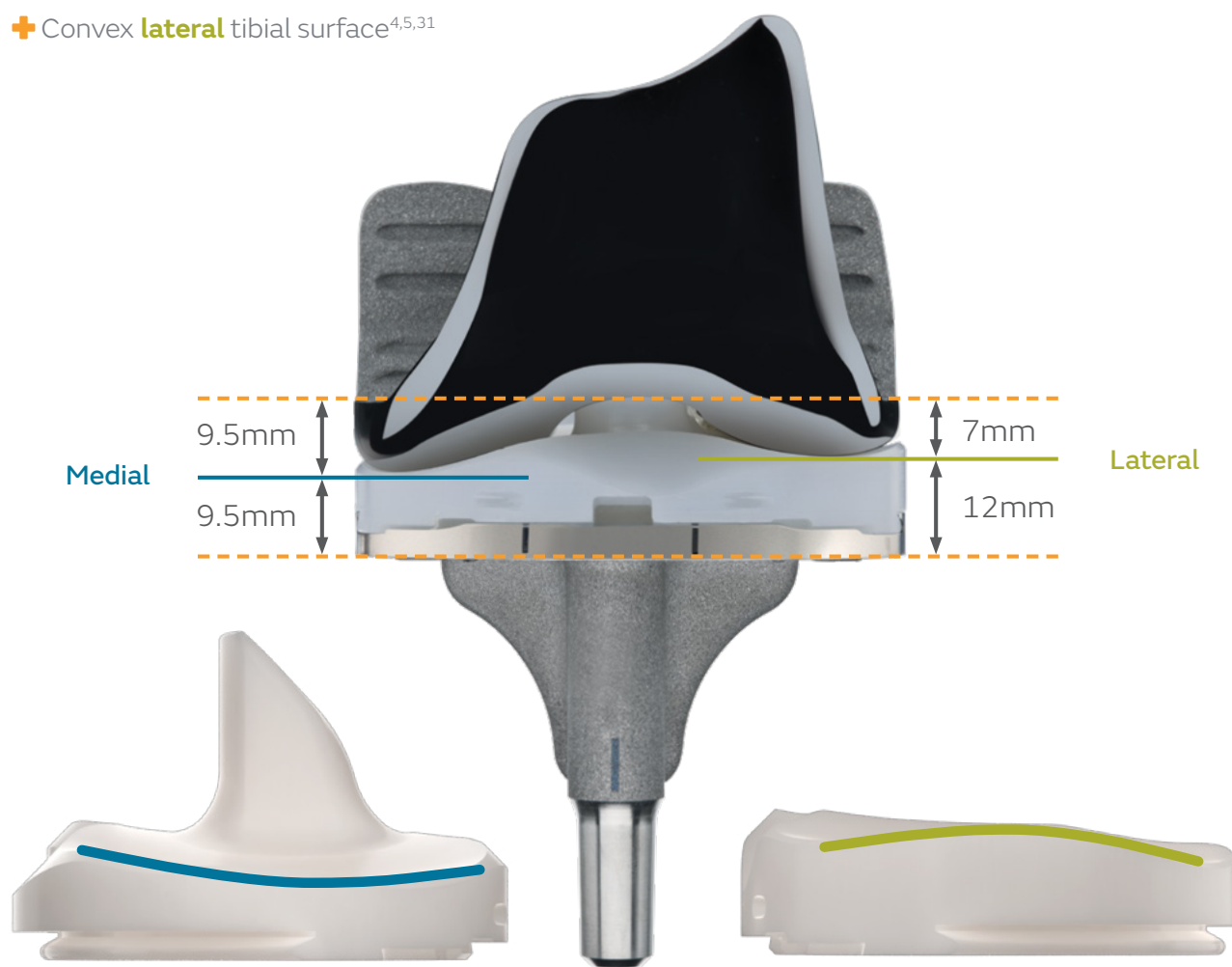

**Medial concavity promotes medial pivot motion patterns<sup>4,5,31,32</sup>**


**Lateral convexity promotes native rollback<sup>4,5,31,32</sup>**

# Normal position<sup>3-5,7,9,31,33-36</sup>


JOURNEY II TKA has also been designed to replicate the mid-line A/P position and 3° varus joint-line found in the normal healthy knee.

The replication of these positions:

- Helps prevent paradoxical motion<sup>4,5,33</sup>
- Promotes muscular efficiency throughout the range of motion<sup>9,33-35</sup>
- Enables natural patellar tracking<sup>+3,7</sup>
- Allows for more normal ligament tension<sup>4,5,36</sup>

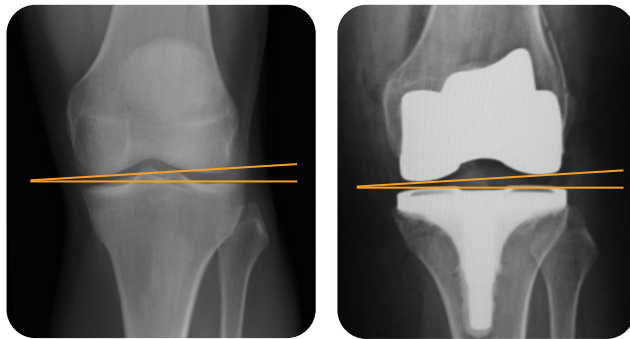

3° Anatomic  
joint-line

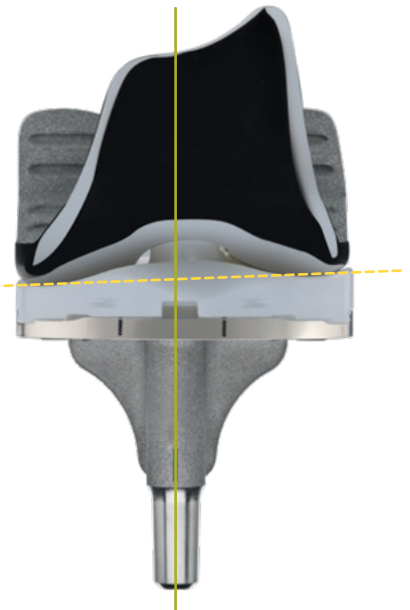

Patella tendon  
angle

No femoral  
overhang




Mid-line  
sulcus  
position




# Normal motion<sup>1-5,8,31,33</sup>

Combining the anatomical shapes and position of the normal knee, JOURNEY II TKA has been able to demonstrate the motion found in the normal knee.


Native knee (left)  
(Grieco)

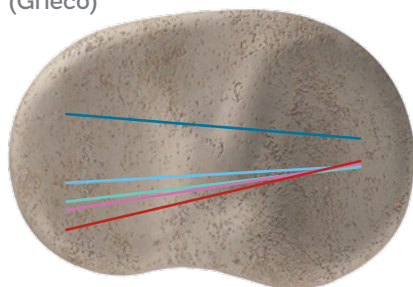

JOURNEY II BCS (left)

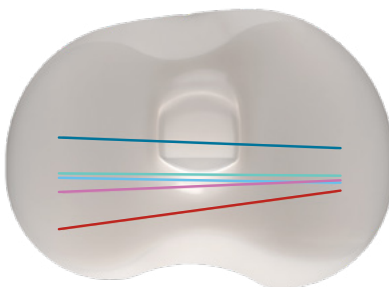



Persona™ MC (left)<sup>20</sup>

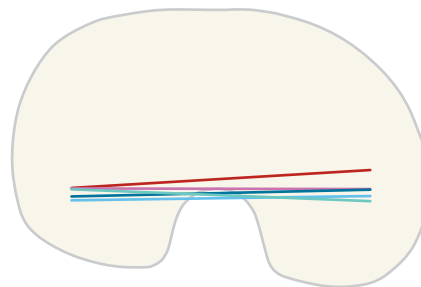

Grieco et al. published that JOURNEY II BCS exhibited:<sup>4</sup>

- Normal-like kinematic patterns and moved as designed under in vivo observation
- Similarities in early and late kinematic patterns with normal asymptomatic knees

Iriuchishima et al. demonstrated that JOURNEY II BCS had no significant difference in rollback ratio or active knee flexion when compared to asymptomatic control knees and Oxford UKA knees.<sup>1</sup>

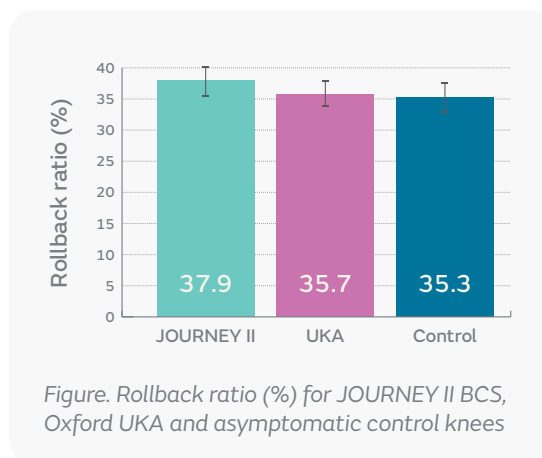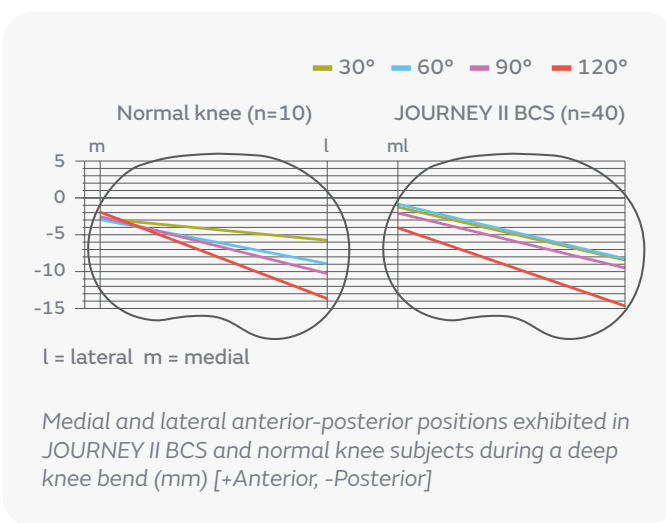

# Rediscover normal: JOURNEY<sup>®</sup> II TKA

The anatomical shape of JOURNEY II TKA is designed to help patients rediscover their normal through a smoother recovery<sup>\*\*6,7</sup> improved function<sup>\*7-11</sup> and higher patient satisfaction.<sup>\*7-9,12</sup>



Smoother recovery



Improved function



Higher patient satisfaction





# Smoother recovery<sup>\*6,7</sup>

## 41%

less likely to be discharged to a skilled nursing facility (p<0.0001)<sup>\*\*+6</sup>

## 35%


more likely to be discharged to home (p=0.0008)<sup>\*\*+6</sup>



Significantly reduced hospital stay (p<0.0001)<sup>\*\*+6</sup>

## Improvements in Mean KSS for JOURNEY II BCS and PS TKA at 6 weeks.

Mean KSS for JOURNEY II BCS and PS TKA<sup>7</sup>

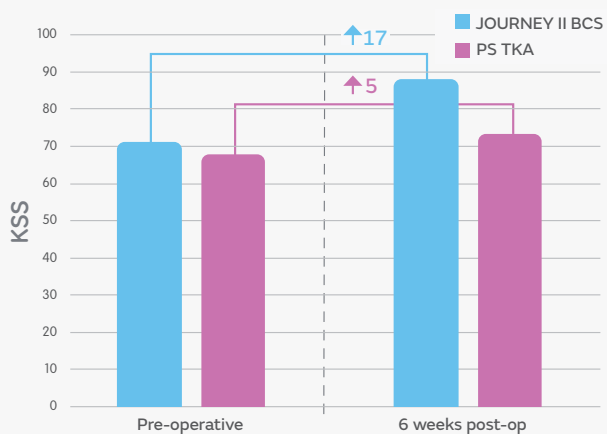

## Improving quadriceps function is important for limiting post-TKA functional deficits<sup>37</sup>


Compared to Attune<sup>®</sup> CR, JOURNEY II CR has shown improved muscle activation and muscle strength in the early recovery period<sup>9</sup>




# Improved function<sup>4-10</sup>

## Gait


JOURNEY II has demonstrated **significant improvements in gait** compared to other TKA designs.<sup>9,34</sup>



Significantly longer step length (p=0.03)<sup>34</sup>




Significantly faster gait speed (p=0.03)<sup>34</sup>




Significantly greater maximum knee extension moment (p=0.04)<sup>34</sup>



Significantly improved knee flexion during walking at 3 months post-TKA (p<0.01)<sup>9</sup>




Significantly improved external rotation during walking at 3 months post-TKA (p<0.01)<sup>9</sup>



Physiological activation timing of investigated muscles at 3 months.<sup>9</sup>

## Flexion

JOURNEY II BCS has shown greater maximal flexion than current TKA designs (Conventional TKA) at 1-year in several studies and a high mean range of motion in a multi-center case study.

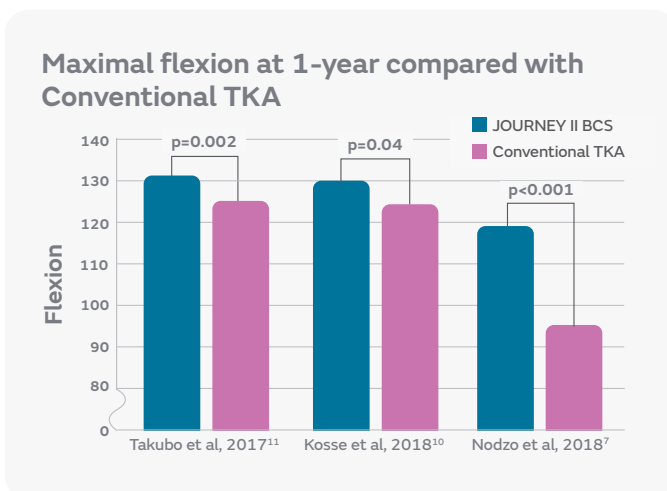




# 131°

mean range of motion at 2 years post-op.<sup>38</sup>

# Higher patient satisfaction<sup>6,7,9,12</sup>

*“Patients reporting their artificial joint as ‘natural’ as opposed to ‘artificial’ are more likely to report higher rates of satisfaction and have higher outcome scores.”<sup>39</sup>*

JOURNEY II has shown significant improvements in Knee Society Scores at 1-year compared to other TKA designs.<sup>7,40</sup>

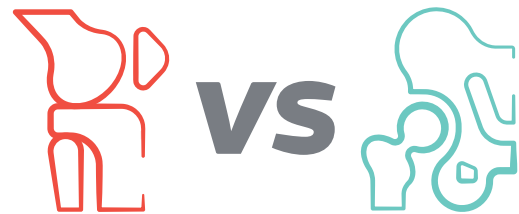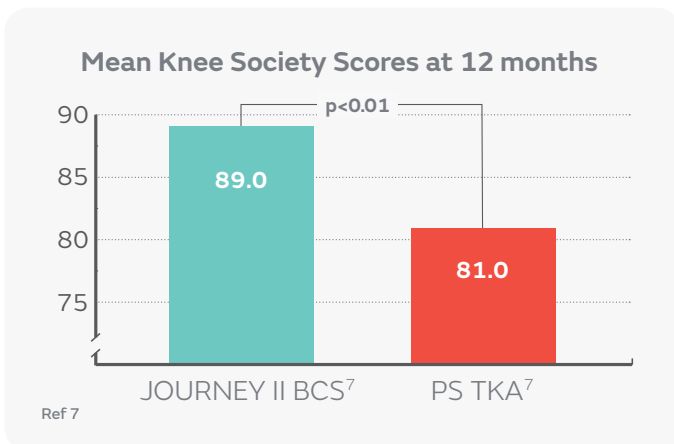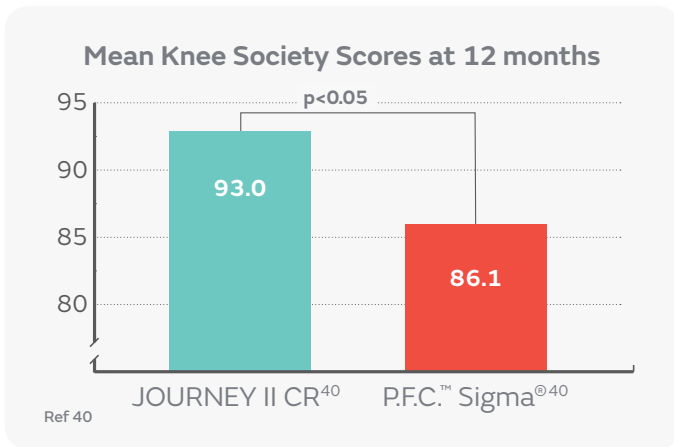

A recent study confirmed that patient satisfaction following TKA remains significantly ( $p = 0.003$ ) lower when compared with THA<sup>39</sup>



Conversely, JOURNEY II has shown **similar outcomes and satisfaction** when compared with clinically similar **THA patients**<sup>41</sup>

## No marked difference in...

- Overall satisfaction at 3 months or 1 year<sup>41</sup>
- Patient quality of life measures at 3 months or 1 year<sup>41</sup>

<sup>‡</sup>Time to return to work, time to return to activities of daily living or time to return to sport

\*Compared to non-JOURNEY II knees  
+Based on BCS evidence

Smith & Nephew Pty Ltd Australia T +61 2 9857 3999 F +61 2 9857 3900	Smith & Nephew Ltd New Zealand T +64 9 820 2840 F +64 9 820 2841	<a href="http://www.smith-nephew.com/australia">www.smith-nephew.com/australia</a> <a href="http://www.smith-nephew.com/new-zealand">www.smith-nephew.com/new-zealand</a> °Trademark of Smith+Nephew. All trademarks acknowledged. 34111-anz V1 REVA 02/25	This material is intended for healthcare professionals. For detailed product information, including indications for use, contraindications, precautions and warnings, please consult the product's applicable Instructions for Use (IFU) prior to use.
---	---	--	--

## References

1. Iriuchishima T, Ryu K. A Comparison of Rollback Ratio between Bicruciate Substituting Total Knee Arthroplasty and Oxford Unicompartmental Knee Arthroplasty. *J Knee Surg.* 2018;31(6):568-572.
2. Murakami K, Hamai S, Okazaki K, et al. Knee kinematics in bi-cruciate stabilized total knee arthroplasty during squatting and stairclimbing activities. *J Orthop.* 2018;15(2):650-654.
3. Carpenter RD, Brillhault J, Majumdar S, Ries MD. Magnetic resonance imaging of in vivo patellofemoral kinematics after total knee arthroplasty. *Knee.* 2009;16(5):332-336.
4. Grieco TF, Sharma A, Dessinger GM, Cates HE, Komistek RD. In Vivo Kinematic Comparison of a Bicruciate Stabilized Total Knee Arthroplasty and the Normal Knee Using Fluoroscopy. *J Arthroplasty.* 2018;33(2):565-571.
5. Smith LA, Nachtrab J, LaCour M, et al. In Vivo Knee Kinematics: How Important Are the Roles of Femoral Geometry and the Cruciate Ligaments? *J Arthroplasty.* 2021;36:1445-1454.
6. Mayman DJ, Patel AR, Carroll KM. Hospital Related Clinical and Economic Outcomes of a Bicruciate Knee System in Total Knee Arthroplasty Patients. Poster presented at: ISPOR Symposium;19-23 May, 2018; Baltimore, Maryland, USA.
7. Nodzo SR, Carroll KM, Mayman DJ. The Bicruciate Substituting Knee Design and Initial Experience. *Tech Orthop.* 2018;33(1):37-41.
8. Murakami K, Hamai S, Okazaki K, et al. In vivo kinematics of gait in posteriorstabilized and bicruciate-stabilized total knee arthroplasties using image-matching techniques. *Int Orthop.* 2018;42(11):2573-2581.
9. Di Benedetto P, Vidi D, Colombo, Buttironi MM, Cainero V, Causero A. Pre-operative and post-operative kinematic analysis in total knee arthroplasty. A pilot study. *Acta Biomed.* 2019;90:91-97.
10. Kosse NM, Heesterbeek PJC, Defoort KC, Wymenga AB, Hellemond G. Minor adaptations in implant design bicruciate-substituted total knee system improve maximal flexion. Poster presented at: 2nd World Arthroplasty Congress;19-21 April, 2018; Rome, Italy.
11. Takubo A, Ryu K, Iriuchishima T, Tokuhashi Y. Comparison of Muscle Recovery Following Bi-cruciate Substituting versus Posterior Stabilized Total Knee Arthroplasty in the Asian Population. *J Knee Surg.* 2017;30(7):725-729.
12. Noble PC, Scuderi GR, Brekke AC, et al. Development of a New Knee Society Scoring System. *Clin Orthop Relat Res.* 2012;470(1):20-32.
13. Noble PC, Gordon MJ, Weiss JM, et al. Does Total Knee Replacement Restore Normal Knee Function? *Clin Orthop Relat Res.* 2005;431:157-165.
14. Scott CEH, Howie CR, MacDonald D, Biant LC. Predicting dissatisfaction following total knee replacement. *J Bone Joint Surg Am.* 2010;92-B(9):1253-1258.
15. Dennis D, Komistek R, Mahfouz M, Walker S, Tucker A. A multicenter analysis of axial femorotibial rotation after total knee arthroplasty. *Clinical Orthopaedics & Related Research.* 2004;428:180-189.
16. Van Onsem S, Verstraete M, Van Eenoo W, Van Der Straeten C, Victor J. Are TKA Kinematics During Closed Kinetic Chain Exercises Associated With Patient-reported Outcomes? A Preliminary Analysis. *Clin Orthop Relat Res.* 2020;478:255-263.
17. Parcells BW, Tria AJ Jr. The Cruciate Ligaments in Total Knee Arthroplasty. *Am J Orthop.* 2016;45:153-160.
18. Collins M, Lavigne M, Girard J, Vendittoli PA. Joint perception after hip or knee replacement surgery. *Orthop Traumatol Surg Res.* 2012;98:275-280.
19. Johal P, Williams A, Wragg P, Hunt D, Gedroyc W. Tibio-femoral movement in the living knee. A study of weight bearing and non-weight bearing knee kinematics using 'interventional' MRI. *J Biomech.* 2005;38(2):269-276.
20. Smith+Nephew 2018. TM-18-064.
21. Smith+Nephew 2018. TM-18-083.
22. Smith+Nephew 2018. TM-18-078.
23. Smith+Nephew 2018. 00225 V3 JOURNEY II TKA Design Rationale 0118.
24. Longo UG, Candela V, Pirato F, et al. Midflexion instability in total knee arthroplasty: a systematic review. *Knee Surg Sports Traumatol Arthrosc.* 2020;10.1007/s00167-020-05909-6. doi:10.1007/s00167-020-05909-6
25. Chambers HG, Sutherland DH. A Practical Guide to Gait Analysis. *J Am Acad Orthop Surg.* 2002;10:222-231.
26. Andriacchi TP, Galante JO, Fermier RW. The Influence of Total Knee-Replacement Design on Walking and Stair-Climbing. *J Bone Joint Surg Am.* 1982;64(9):1328-1335.
27. Dorr LD, Ochsner JL, Gronley J, Perry J. Functional Comparison of Posterior Cruciate-Retained versus Cruciate-Sacrificed Total Knee Arthroplasty. *Clin Orthop Relat Res.* 1988;236:36-43.
28. Kramers-de Quervain IA, Stussi E, Muller R, et al. Quantitative Gait Analysis After Bilateral Total Knee Arthroplasty With Two Different Systems Within Each Subject. *J Arthroplasty.* 1997;12(2):168-179.
29. Saari T, Tranberg R, Zügner R, Uvehammer J, Kärrholm J. Changed gait pattern in patients with total knee arthroplasty but minimal influence of tibial insert design. *Acta Orthop.* 2005;76(2):253-260.
30. Lester DK, Shantharam R, Zhang K. Dynamic electromyography after cruciate-retaining total knee arthroplasty revealed a threefold quadriceps demand compared with the contralateral normal knee. *J Arthroplasty.* 2013;28(4):557-562.
31. Kaneko T, Kono N, Mochizuki Y, Hada M, Toyoda S, Mushi Y. Bicruciate substituting total knee arthroplasty improved medio-lateral instability in mid-flexion range. *J Orthop.* 2017;14(1):201-206.
32. Brillhault J, Ries MD. Measuring patellar height using the lateral active flexion radiograph: Effect of total knee implant design. *Knee.* 2010;17(2):148-151.
33. Catani F, Ensini A, Belvedere C, et al. In vivo kinematics and kinetics of a bi-cruciate substituting total knee arthroplasty: a combined fluoroscopic and gait analysis study. *J Orthop Res.* 2009;27(12):1569-1575.
34. Hyodo K, Kanamori A, Kadone H, Takahashi T, Kajiwara M, Yamazaki M. Gait Analysis Comparing Kinematic, Kinetic, and Muscle Activation Data of Modern and Conventional Total Knee Arthroplasty. *Arthroplast Today.* 2020;6(3):338-342.
35. Hada M, Mizu-Uchi H, Okazaki K, et al. Bi-cruciate stabilized total knee arthroplasty can reduce the risk of knee instability associated with posterior tibial slope. *Knee Surg Sports Traumatol Arthrosc.* 2018;26(6):1709-1716.
36. Smith+Nephew 2012. JRN2 KneeSim Analysis Memo.
37. Thomas AC, Stevens-Lapsley JE. Importance of attenuating quadriceps activation deficits after total knee arthroplasty. *Exerc Sport Sci Rev.* 2012;40:95-101.
38. Harris AI, Luo TD, Lang JE, Kopjar B. Short-term safety and effectiveness of a second-generation motion-guided total knee system. *Arthroplast Today.* 2018;4:240-243.
39. Varacallo M, Chakravarty R, Denehy K, Star A. Joint perception and patient perceived satisfaction after total hip and knee arthroplasty in the American population. *J Orthop.* 2018;15:495-499.
40. Lutes W, Fitch D. Comparison of functional outcomes following total knee arthroplasty with a conventional implant design or one designed to mimic natural knee kinematics. Presented at: 39th SICOT Orthopaedic World Congress; October 10-13, 2018; Montréal, Canada.
41. Snyder MA, Sympon A, Gregg J, Levit A. A comparison of patient reported outcomes between total knee arthroplasty patients receiving the Journey II bi-cruciate stabilizing knee system and total hip arthroplasty patients. *Ortop Travmatol Protez.* 2018;3:5-10.