


+ OneNPWT clinical decision tree for open wounds

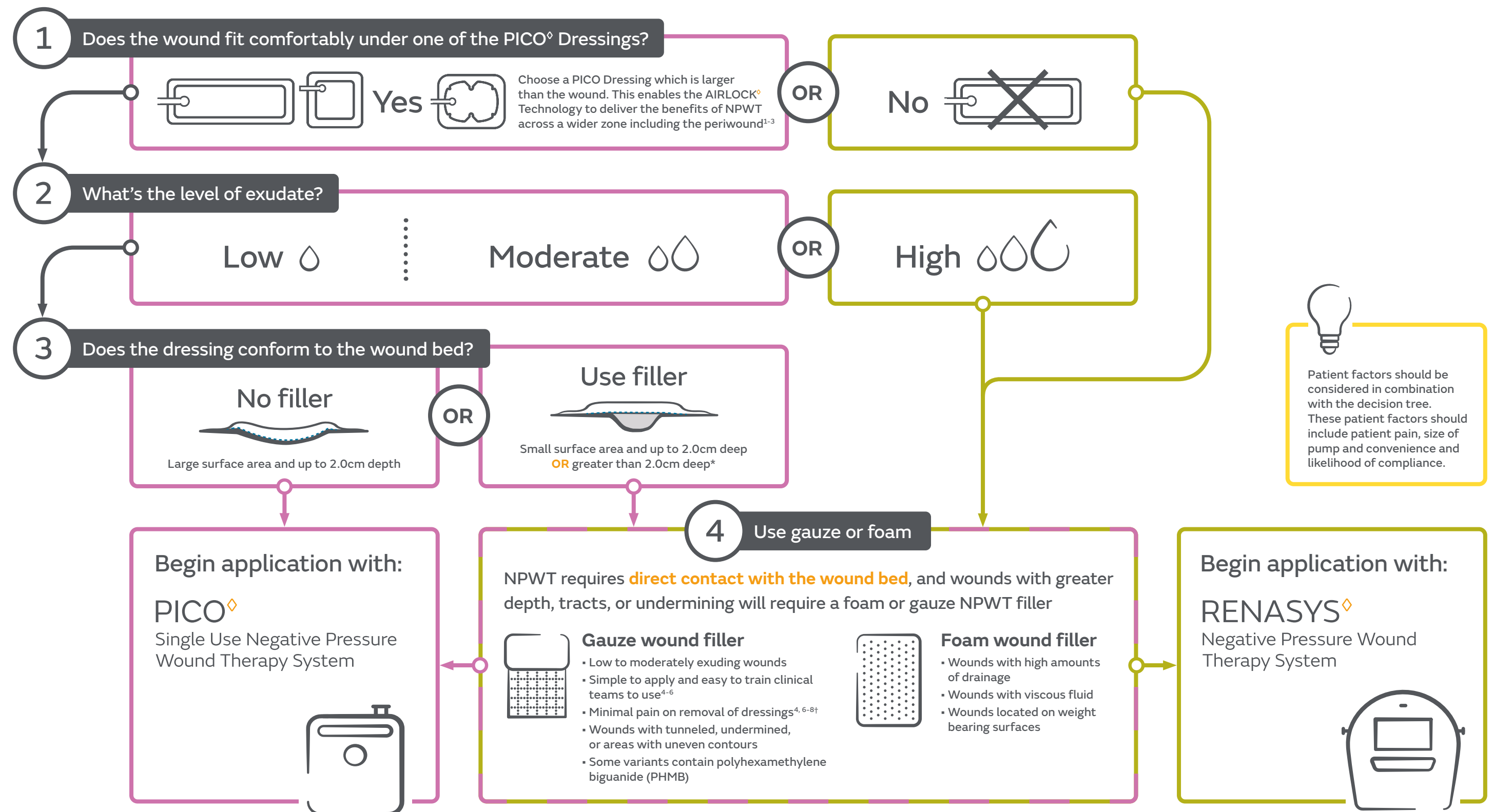

Smith+Nephew does not provide medical advice. The information presented is not, and is not intended to serve as, medical advice. It is the responsibility of healthcare professionals to determine and utilize the appropriate products and techniques according to their own clinical judgment for each of their patients. The information presented may not be appropriate for all jurisdictions. For detailed product information, including indications for use, contraindications, precautions and warnings, please consult the product's applicable Instructions for Use (IFU) prior to use.

* Wounds must not contain exposed arteries, veins, nerves or organs. † p=0.046; n=31; Compared to black foam in acute post traumatic wounds. **Reference:** **1.** Smith & Nephew 2019.PICO Biomechanical Study. Internal Report. DS/19/211/R. **2.** Casey C. Consistent delivery of therapeutic negative pressure levels by a single use negative pressure wound therapy system (sNPWT)* in a wound model. Paper presented at: EWMA; 2019; Gothenburg, Sweden. **3.** Smith+Nephew 2021.PICO™ Pressure Mapping Study. Internal Report. DS/19/211/R - Part B. **4.** Hurd T, Chadwick P, Cote J, Cockwill J, Mole T, Smith J. Impact of gauze-based NPWT on the patient and nursing experience in the treatment of challenging wounds. International Wound Journal. 2010;7(6):448-455. **5.** Fracalvieri M, Scalise A, Ruka E, et al. Negative pressure wound therapy using gauze and foam: Histological, immunohistochemical, and ultrasonography morphological analysis of granulation and scar tissues – Second phase of a clinical study. In. European Journal of Plastic Surgery. Vol 37 2014:411-416. **6.** Johnson S. VISTA® – A new option in Negative Pressure Therapy. Journal of Wound Technology. 2008;1:30-31. **7.** Fracalvieri M, Ruka E, Bocchiotti M, Zingarelli E, Bruschi S. Patient's pain feedback using negative pressure wound therapy with foam and gauze. International wound journal. 2011;8(5):492-499. **8.** Smith+Nephew 2009. A prospective, open labelled, multicentre evaluation of the use of VISTA in the management of chronic and surgical wounds and A prospective, open labelled evaluation of the use of EZCare in the management of chronic and acute wounds. Internal Report. SR/CIME/010/012. [◇]Trademark of Smith+Nephew. All Trademarks acknowledged. ©July 2023 Smith+Nephew. AWM-AWD-38307 | GMC1146d | RoW