


GoFlo Pump

Simplicity without compromise

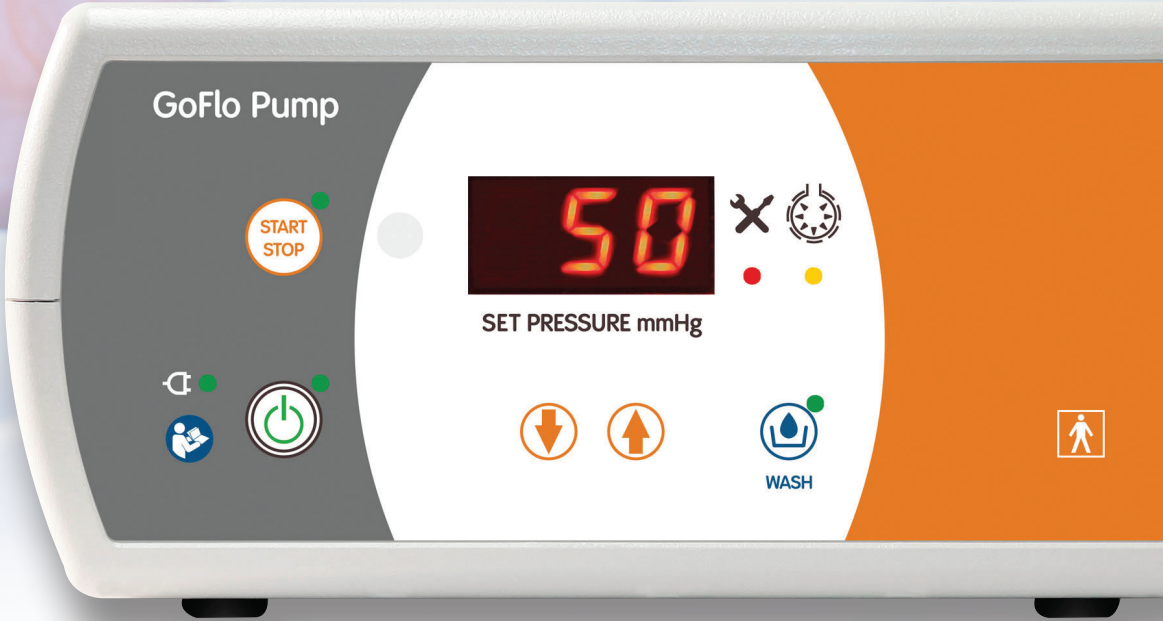

Effective visualization for a broad range of procedures

DESIGNED TO SIMPLY GET THE JOB DONE...

Consistently

- Flow regulation is optimized through automatic cannula detection
- Designed to maintain distention throughout procedure
- Continuous self-calibration ensures settings are reproducible

Effectively

- Broad surgical capability for all joint procedures
- Optimal flow performance using any scope cannula

Intuitively

- Simple user interface and wireless remote
- Quick and uncomplicated tube set connection
- Small footprint, either pole mounted or stacked
- Automatic cannula detection executes in first 20 seconds of run cycle

FEATURES & SPECIFICATIONS

Automatic self-test at startup

- Ensures the pump is running correctly and is optimized for specific instruments or procedures
- 20 second initial calibration cycle to adjust performance based on cannula specifications
 - Recommended to perform outside the joint space but can be conducted within the joint space

Specifications

- Max flow: 1.5 l/min
- Pressure range: 5-150mmHg
- RFID tube set recognition – prevents the use of counterfeit tubing or reprocessing
- “Wash” function – increase in pressure +50% for 20 seconds

Ordering Information

GoFlo Pump Reference #	Description
72204968	GoFlo Pump
T0449-01	Sterile Tubeset 10/box
72204972	GoFlo Remote Control
72205012	Remote Control Sterile Sleeve 50/box

Manufactured by:
W.O.M. World of Medicine GmbH

Distributed by:
Smith & Nephew, Inc. www.smith-nephew.com
150 Minuteman Road T +1 978 749 1000
Andover, MA 01810 USA US Customer Service: +1 800 343 5717



Setup is easy as:



1



● STEP 1: Power On

2



● STEP 2: Load Tubing

3



● STEP 3: Start and Go