

+ Position with ease

The beach chair positioner that
sits well with patients, surgeons,
anesthesiologists and OR staff

Smith+Nephew



T-MAX
Beach Chair



Access and versatility Made easy

Our T-MAX Beach Chair positioner provides access, stability and versatility for arthroscopic and open shoulder procedures in the beach chair position.¹ It places the surgeon, anesthesiologist and surgical staff in an optimal position to position patients.

The chair's adjustable legs allow the T-MAX to attach to many commonly used OR tables and securely positions patients of up to 350 lbs, with a maximum height of 6'5" (1.96m).² Our T-MAX securely holds the patient while positioning and re-positioning between 0 and 90 degrees and the built-in lateral supports provide sturdy counter traction when distracting the operative arm. With a dual ball and socket head rest, adjustable lateral supports, sliding shoulder pads and more, the T-MAX is designed to position the patient with ease and enables the surgeon unrestricted access to the surgical site.¹

+ Hydraulic Assist Mechanism



KEY FEATURES

+ Intubation tree (static)

Intubation tree is used to hold intubation tubing and can simply be rotated out of the way if not required. An additional head positioning option providing anesthesiologists with a broad, flat surface for intubations

+ Head rest

Double ball and socket head rest offers infinite positions of tilt and rotation, maintaining the patient's head and airway in a static position throughout the procedure

With wide range of tilt and rotation possibilities, the anesthesiologist can choose their desired head position



Head rest, extended



Head rest, retracted



Head rest angled for left shoulder access



Head rest angled for right shoulder access

+ Shoulder pad

The shoulder pad is designed to provide stability. Sliding shoulder pad and posterior placement can allow for unrestricted access to the surgical site¹



Sliding Shoulder pad positioned to provide support for shoulder arthroplasty



Sliding Shoulder pad positioned to provide full exposure to the medial boarder of the scapula for shoulder arthroscopy

+ Non-operative arm board (static)

Provides support and adjustability for the non-operative arm, eliminating the need to attach an arm board elsewhere to the OR table

+ Lateral Adjustment

Lateral adjustment allows the T-MAX to slide laterally on the frame to either side of the OR table, to position patients close to the surgeon

+ Spider Rail

The only beach chair designed to work with the SPIDER Limb Positioner. The T-MAX system contains a SPIDER Rail which allows the SPIDER to be placed on the back of the T-MAX offering the surgeon additional space in the surgical field

+ Lateral support pads

The Lateral Support Pads are designed to ensure patient stability while providing counter traction when distracting the operative arm

+ Hydraulic assist mechanism

Hydraulic assist mechanism assists in raising and lowering patient to any angle between 0° and 90°



Ordering information

T-MAX Beach Chair Positioner

Reference #	Description
7210551	T-MAX Beach Chair Positioner Includes
	Chair and Pad Assembly
	Head Positioning Assembly
	Table Leg Assembly
	Knee Bolster
	Lateral support (2 Units)

Disposable Accessories

7210559	Disposable Face Masks (case of 8)
---------	-----------------------------------


T-MAX Table Rail Clamps and Leg Kits

Reference #	Description
7210554	North American Rail Clamps (2 Units)
7211099	European Square Rails Clamps (2 Units)
7211096	JE Adaptor Leg Kit
7211097	JE Rail Clamps (2 Units)
7211098	DE Adapter Leg Kit
7211100	Denyer Rail Clamps (2 Units)

T-MAX Accessories

Reference #	Description
7210552	Adjustable Arm Board
7210553	T-MAX Cart
7210560	Knee Bolster
7210563	Drape Bar
7210564	Lateral Support Knobs (2 Units)
7210555	Intubation Pad
7210556	Gel Face Mask
7210557	Extended Lateral Support (2 Units)
7210558	Mask Straps (pair)
7210568	Table Stabilization Straps
72203415	Non-Operative Arm Board
72203318	Adjustable Lateral Support

Highlight T-MAX Accessories

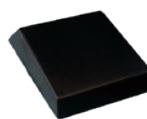

7210560 – Knee Bolster
The knee bolster helps prevent the patient from sliding down the OR table.



7210563 – Drape Bar



7210555 – Intubation Pad
The intubation pad (optional) provides anesthesiologists with a broad, flat surface for ease-of-use with intubation.



7210559 – Disposable Face Masks (case of 8)



72203415 – Non-operative Arm Board



72203318 – Adjustable Lateral Support

Smith & Nephew Pty Ltd
Australia
T +61 2 9857 3999
F +61 2 9857 3900

Smith & Nephew Ltd
New Zealand
T +64 9 820 2840
F +64 9 820 2841

www.smith-nephew.com/australia
www.smith-nephew.com/new-zealand
°Trademark of Smith+Nephew. All trademarks acknowledged. 02942-anz V3 05/21

For detailed product information, including indications for use, contraindications, precautions and warnings, please consult the product's applicable Instructions for Use (IFU) prior to use.

References

1. Verhaegen F, Debeer P. *Arthroscopic Removal of Rotator Cuff Calcifications: Operative Technique*. JBJS essential surgical techniques. 2016;6(4):e38-e38. 2. Data on file: Table test calculations document updated 8/19/2014.