



+ High exudate?

Let's  
deal with it

With the no-compromise fluid management of DURAMAX<sup>◇</sup> N  
Non-Adhesive Super Absorbent Dressing.<sup>1</sup> **Just when you need it.**

**Smith+Nephew**



DURAMAX<sup>◇</sup> N  
Non-Adhesive  
Super Absorbent Dressing

# Performance, protection, and maximum absorption...



## MAXIMUM ABSORPTION

DURAMAX<sup>o</sup> N Non-Adhesive Super Absorbent Dressing manages excess exudate:



### ABSORPTION

The dressing **absorbs and locks in** viscous exudate.<sup>1</sup>



### CONFORMABLE

Can help to **stay in close contact** with skin.




### COMPRESSION

The dressing can absorb and retain fluid **while under standard compression.**<sup>1\*</sup>



## DEPENDABLE PROTECTION

DURAMAX N Dressing absorbs a large amount of exudate, and is comfortable for patients:

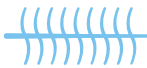

### BREATHABLE<sup>1</sup>

The dressing is breathable and backing layer acts as a **fluid repellent** on outer surface.<sup>2</sup>



### BACTERIA

Once within the dressing, **over 99% of bacteria are locked within the dressing** and away from the wound.<sup>3†</sup>



### RETENTION CAPACITY

The dressing **minimises risk of maceration and irritation of the wound bed** and periwound tissue due to exudate.<sup>1</sup>

\*40mmHg. †Tested in vitro on *Pseudomonas aeruginosa* and Methicillin-resistant *Staphylococcus aureus*.



## OPTIMAL PERFORMANCE

DURAMAX<sup>◊</sup> N Non-Adhesive Super Absorbent Dressing is for acute and chronic wounds with **moderate to heavy exudate**. Typical wounds are:



DIABETIC  
FOOT ULCERS



PRESSURE  
ULCERS



VENOUS LEG &  
ARTERIAL ULCERS



FIRST & SECOND  
DEGREE BURNS



POST-OP  
WOUNDS



TRAUMATIC  
WOUNDS

DURAMAX N Dressing features a **5-layer** structure:

### SMS non-woven backing layer

- Breathable and water repellent<sup>1,2</sup>

### Superabsorbent core layer

- Absorb and retain exudate<sup>1</sup>

### Hydrophilic polypropylene (PP) non-woven layer

- Absorb exudate, transfer fast and evenly<sup>1</sup>

### Hydrophilic polyethylene (PE)/ polyethylene terephthalate (PET) contact layer

- Comfortable wound contact layer, absorb and transfer exudate<sup>4</sup>



... when you need it


# DURAMAX<sup>◊</sup> N Non-Adhesive Super Absorbent Dressing is for moderate to heavy exudate

For moderate to heavy exudate




## DURAMAX N Non-Adhesive Super Absorbent Dressing — fluid retention

DURAMAX<sup>◊</sup> N Dressing has demonstrated significantly\* superior fluid retention versus leading competitor.<sup>5</sup>

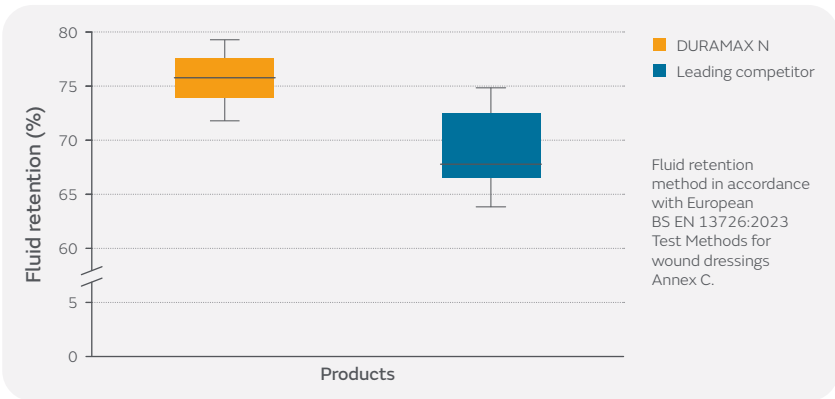

## DURAMAX N Non-Adhesive Super Absorbent Dressing — absorption during compression

DURAMAX<sup>◊</sup> N Dressing has demonstrated significantly† superior absorption during compression versus other super absorbent dressings.<sup>6</sup>

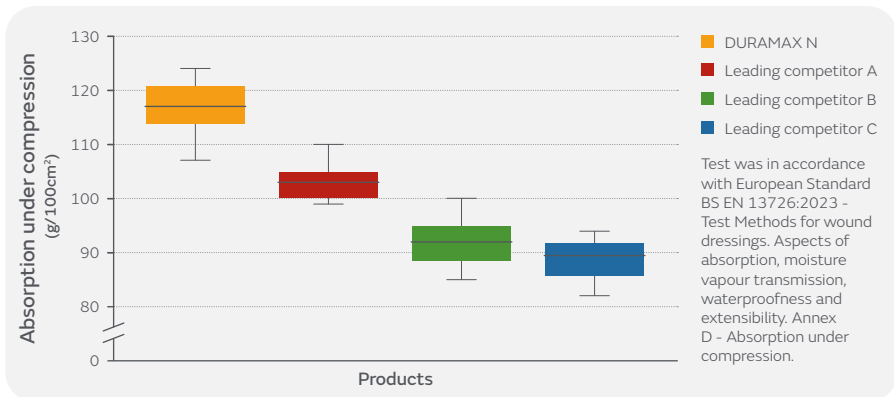

\*p<0.01. †p<0.0167.

## DURAMAX<sup>o</sup> N Dressing is easy to apply

1



Select a dressing size that is **big enough to overlap the wound edges** by at least 2-3cm.

2



**Gently apply the dressing** directly on the wound site with the white hydrophilic side of dressing facing the wound surface.

3



The dressing needs to be **secured by an appropriate secondary dressing**, such as adhesive tape an elastic bandage or a pressure bandage.



## Ordering DURAMAX<sup>◇</sup> N Dressing

Code	Size	Qty
66807730	10cm x 10cm	10
66807731	15cm x 20cm	10
66807732	20cm x 20cm	10
66807733	20cm x 25cm	10
66807734	20cm x 40cm	10



### DURAMAX<sup>◇</sup> N

Non-Adhesive  
Super Absorbent Dressing

This material is intended for healthcare professionals. For detailed product information, including indications for use, contraindications, precautions and warnings, please consult the product's applicable **Instructions for Use (IFU)** prior to use.

**References:** **1.** Test Report NO:JCSQ2025021910 from Winner Medical Lab. **2.** Test Report NO: JY2-241102-02 from Winner Medical Lab. **3.** Test Report NO: JY2025040701 from Winner Medical Lab. **4.** Test report NO:JY20250410-01 from Winner Medical Lab. **5.** Smith+Nephew 2025. Internal report CSD.AWM.25.038. **6.** Smith+Nephew 2025. Internal report CSD.AWM.25.037.

#### Advanced Wound Management

Smith & Nephew Pty Ltd (Australia) +61 2 9857 3999 [smith-nephew.com/en-au](http://smith-nephew.com/en-au).

Smith & Nephew Ltd (New Zealand) +64 9 820 2840 [smith-nephew.com/en-nz](http://smith-nephew.com/en-nz).

◇Trademark of Smith+Nephew. All trademarks acknowledged. © April 2026

Smith+Nephew. 47727-anz V1 04/26.

[www.smith-nephew.com](http://www.smith-nephew.com)