Simple. Flexible. Gentle.

The RENASYS Adhesive Gel Patch is a sterile, single use patch intended for fixation of drainage tubing and to help improve seals in support of existing Smith+Nephew NPWT dressing kits. Compatible with PICO◊ Single Use Negative Pressure Wound Therapy and RENASYS◊ Negative Pressure Wound Therapy.

Simple – Helps busy clinicians easily address leaks in complex areas and difficult-to-seal wounds

Smith+Nephew

RENASYS◊
Adhesive Gel Patch

Product number	Description	Size	Patches/box	Boxes/case
66801082	RENASYS Adhesive Gel Patch	7cm x 10cm	10	5



Application techniques: creating a seal

Challenging anatomical areas



Cut the Gel Patch into strips in a direction with backing removal ends accessible.




Remove the backing on one side only and apply to skin with gentle pressure. Ensure the fold is addressed first.




Remove the backing. Apply Transparent Film over the foam or gauze interface and the Gel Patch to create a seal and finish the dressing.

Around wound margins

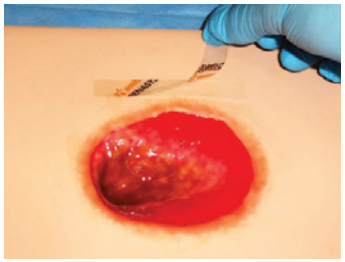

Cut Gel Patch into several strips. Remove the backing on one side only and apply to skin with gentle pressure. Remove the backing.

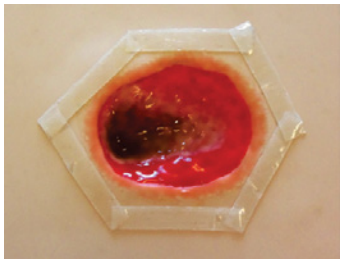

Continue to apply Gel Strips around wound margins ensuring that the strips overlap to create a good seal. Continue with normal dressing application.

Challenging areas of the foot



Cut the Gel Patch into strips in a direction with backing removal ends accessible.





Remove the backing on one side only and apply to skin with gentle pressure. Remove the backing. Apply Transparent Film over the foam or gauze interface and the Gel Patch to create a seal and finish the dressing.

Advanced Wound Management
Smith & Nephew Inc.
Fort Worth, TX 76109 USA

Customer Care Center

1-800-876-1261
F 727-392-6914

www.smith-nephew.com

©2025 Smith+Nephew, Inc.
◊All Trademarks acknowledged.

RNBE6-46538-0225